

Estradiol - USP Method (Synergi 4 μ m Fusion-RP)**Column:** SynergiTM 4 μ m Fusion-RP 80 Å, LC Column 150 x 4.6 mm, Ea**Dimensions:** 150 x 4.6 mm ID**Order No:** 00F-4424-E0**Elution Type:** Isocratic**Eluent A:** Water**Eluent B:** Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	45	55

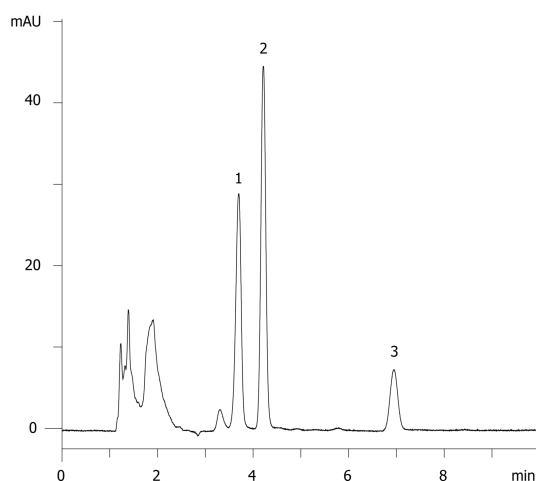
Flow Rate: 1 mL/min**Col. Temp.:** 25 °C**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 205 nm (25 °C)**Analyst Note:** SecurityGuardTM Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Fusion-RP 4 x 3.0mm, 10/Pk Part No.: AJ0-7557

- Holder Part No.: KJ0-4282



15663

**ANALYTES:****1** Ethyl paraben

Retention Time: 3.697 min

2 Estradiol

Retention Time: 4.213 min

3 Estrone

Retention Time: 6.946 min

