

## Ginseng (Gemini 5µm C18)

**Column:** Gemini® 5 µm C18 110 Å, LC Column 150 x 4.6 mm, Ea

**Dimensions:** 150 x 4.6 mm ID

**Order No:** 00F-4435-E0

**Elution Type:** Gradient

**Eluent A:** Water

**Eluent B:** Acetonitrile

| Gradient Profile: | Step No. | Time (min) | Pct A | Pct B |
|-------------------|----------|------------|-------|-------|
|                   | 1        | 0          | 85    | 15    |
|                   | 2        | 60         | 20    | 80    |

**Flow Rate:** 1.5 mL/min

**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 205 nm (ambient)

**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Gemini C18 4 x 3.0mm, 10/Pk Part No.: AJ0-7597


- Holder Part No.: KJ0-4282



Products used in this application:



15693

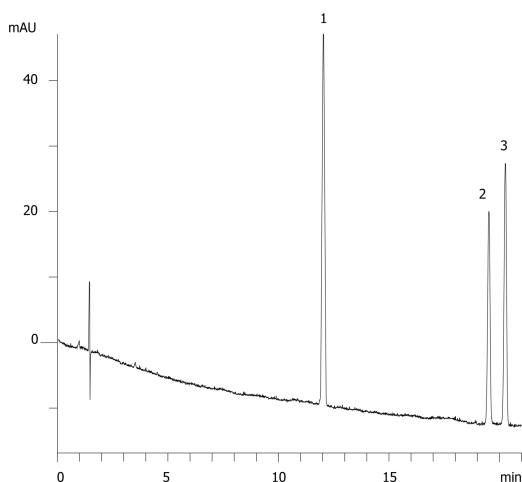

### ANALYTES:

- 1 Ginsenoside-Re  
Retention Time: 12.02 min
- 2 Ginsenoside-Rb1  
Retention Time: 19.52 min
- 3 Ginsenoside-Rc  
Retention Time: 20.27 min

