

Organics on Luna 3u C18(2) Capillary LC

Column: Luna® 3 µm C18(2) 100 Å, LC Column 50 x 0.3 mm, Capillary, Ea

Dimensions: 50 x 0.3 mm ID

Order No: 00B-4251-AC

Elution Type: Gradient

Eluent A: 10mM Ammonium Acetate

Eluent B: Acetonitrile

Gradient	Step No.	Time (min)	Pct A	Pct B
----------	----------	------------	-------	-------

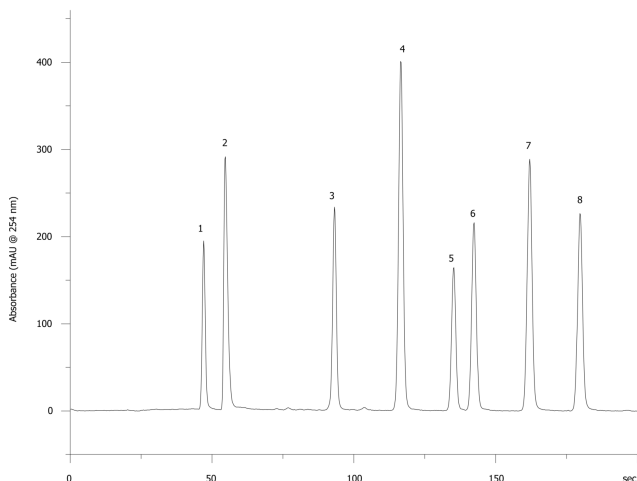
Profile:	1	0	98	2
	2	4	2	98

Flow Rate: 10 $\mu\text{L}/\text{min}$

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

15720

**ANALYTES:**

- 1 Theophylline
- 2 5-Phenyl-1H-Tetrazole
- 3 Colchicine
- 4 Acetophenone
- 5 Indole
- 6 Propiophenone
- 7 Butyrophenone
- 8 Valerophenone



Products used in this application:

