

Cytosine family using Synergi Polar-RP (acetate)**Column:** SynergiTM 4 μ m Polar-RP 80 Å, LC Column 150 x 4.6 mm, Ea**Dimensions:** 150 x 4.6 mm ID**Order No:** 00F-4336-E0**Elution Type:** Gradient**Eluent A:** 20mM Ammonium acetate, pH 4.5**Eluent B:** Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	97	3
	2	10	70	30

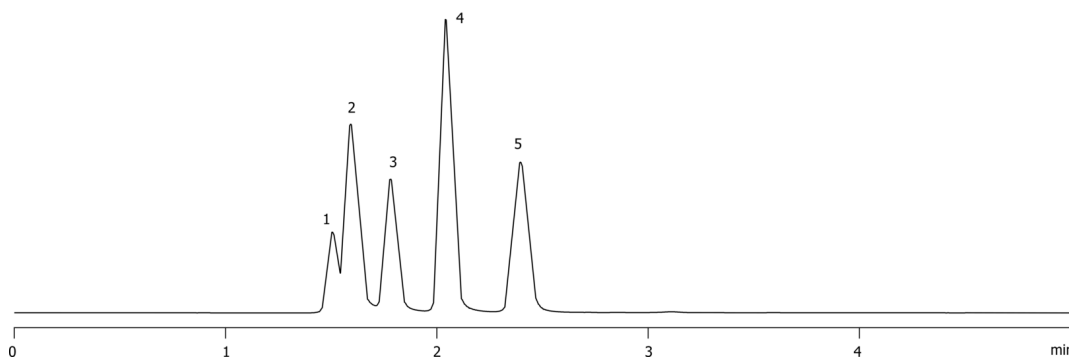
Flow Rate: 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 260 nm (ambient)**Analyst Note:** SecurityGuardTM Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, POLAR-RP 4 x 3.0mm, 10/Pk Part No.: AJ0-6076

- Holder Part No.: KJ0-4282



15762

**ANALYTES:**

- 1 CTP
- 2 CDP
- 3 CMP
- 4 Cytosine
- 5 Cytidine

