

Uracil family using Synergi Polar-RP (acetate)

Column: SynergiTM 4 μ m Polar-RP 80 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4336-E0

Elution Type: Gradient

Eluent A: 20mM Ammonium acetate, pH 4.5

Eluent B: Acetonitrile

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	97	3
	2	10	80	20

Flow Rate: 1 mL/min

Col. Temp.: ambient


Detection: UV-Vis Abs.-Variable Wave.(UV) @ 260 nm (ambient)

Analyst Note: SecurityGuardTM Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, POLAR-RP 4 x 3.0mm, 10/Pk Part No.: AJ0-6076

- Holder Part No.: KJ0-4282

15764

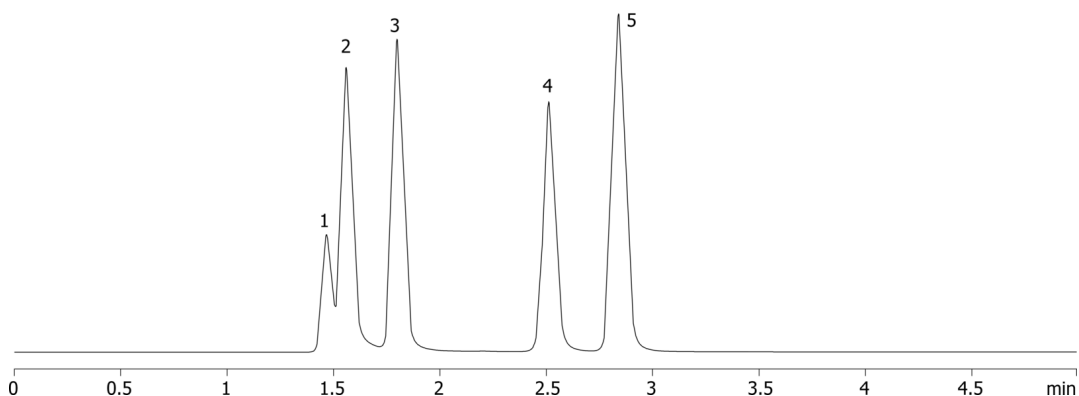

ANALYTES:

- 1** UTP
- 2** UDP
- 3** UMP
- 4** Uracil
- 5** Uridine

