

Deoxynucleosides using Synergi Fusion-RP

Column: SynergiTM 4 μ m Fusion-RP 80 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4424-E0

Elution Type: Gradient

Eluent A: 20mM Ammonium acetate, pH 4.5

Eluent B: Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	10	60	40

Flow Rate: 1 mL/min

Col. Temp.: ambient


Detection: UV-Vis Abs.-Variable Wave.(UV) @ 260 nm (ambient)

Analyst Note: SecurityGuardTM Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Fusion-RP 4 x 3.0mm, 10/Pk Part No.: AJ0-7557

- Holder Part No.: KJ0-4282

15766

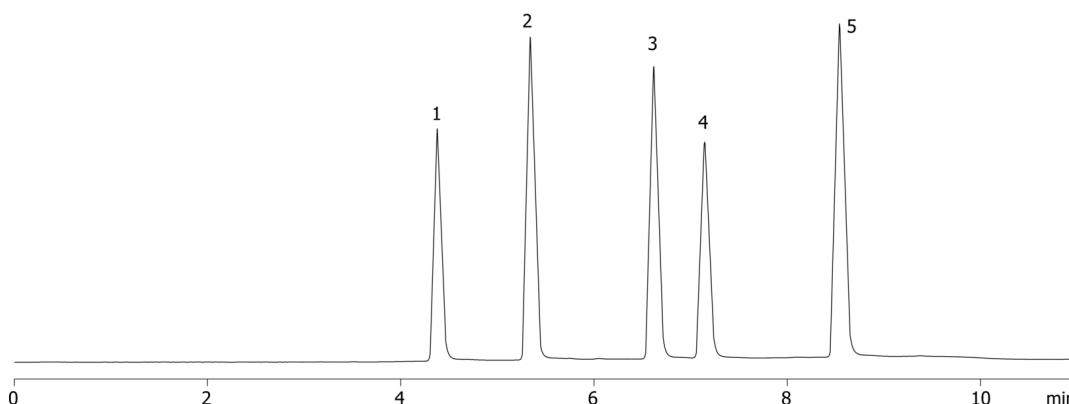

ANALYTES:

- 1 2'-Deoxycytidine
- 2 2'-Deoxyuridine
- 3 2'-Deoxyguanosine
- 4 2'-Deoxythymine
- 5 2'-Deoxyadenine

