

Nucleosides using Synergi Fusion-RP**Column:** SynergiTM 4 μ m Fusion-RP 80 Å, LC Column 150 x 4.6 mm, Ea**Dimensions:** 150 x 4.6 mm ID**Order No:** 00F-4424-E0**Elution Type:** Gradient**Eluent A:** 20mM Ammonium acetate pH 4.5**Eluent B:** Acetonitrile

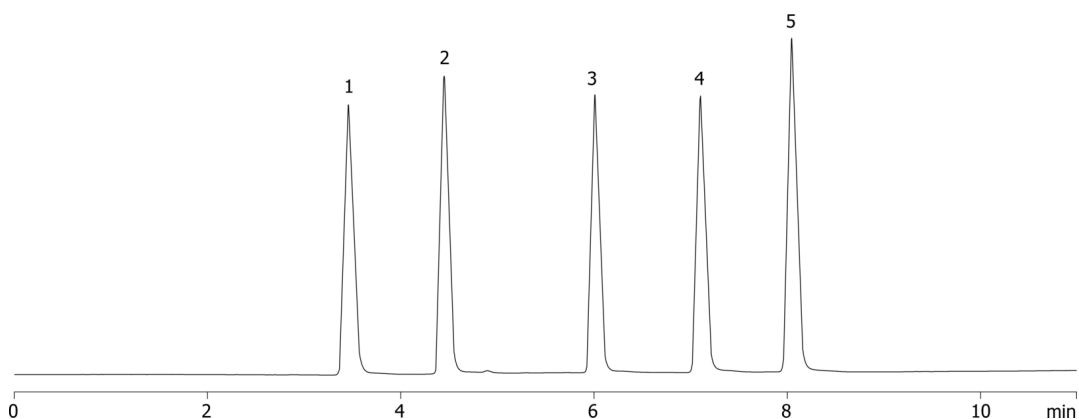
Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	10	60	40

Flow Rate: 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 260 nm (ambient)**Analyst Note:** SecurityGuardTM Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Fusion-RP 4 x 3.0mm, 10/Pk Part No.: AJ0-7557

- Holder Part No.: KJ0-4282

15768

**ANALYTES:**

- 1 Cytidine
- 2 Uridine
- 3 Guanosine
- 4 Thymine
- 5 Adenine

