

Sialic acid mix using Rezex ROA

Column: Rezex™ ROA-Organic Acid H+ (8%), LC Column 300 x 7.8 mm, Ea

Dimensions: 300 x 7.8 mm ID

Order No: 00H-0138-K0

Elution Type: Isocratic

Eluent A: 0.08% Acetic acid in Water

Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100

Flow Rate: 0.5 mL/min

Col. Temp.: 80 °C

Detection: ELSD @ (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Carbo-H 4 x 3.0mm, 10/Pk Part No.: AJ0-4490

- Holder Part No.: KJ0-4282


15769

REZEX™
HPLC COLUMNS for Carbohydrates, Oligosaccharides, and Organic Acid Analysis



Products used in this application:

REZEX™

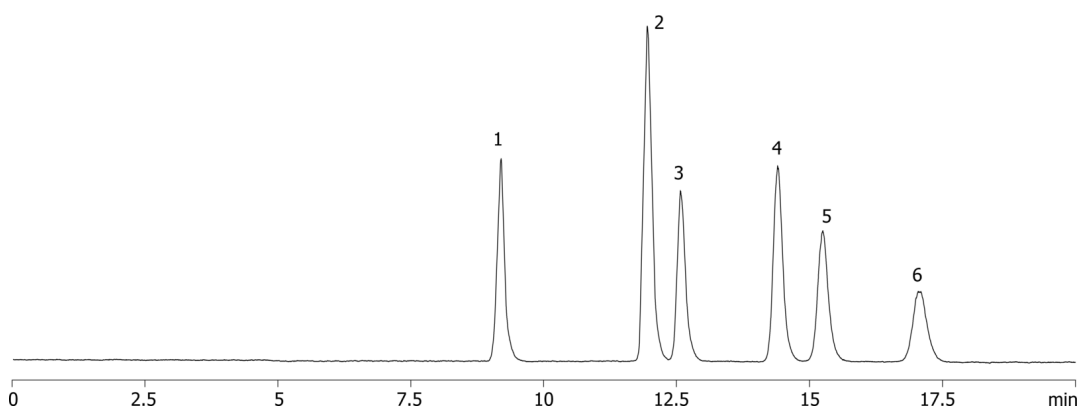

ANALYTES:

- 1 Sialic acid
- 2 Glucose
- 3 Mannose
- 4 Fucose
- 5 N-Acetyl-D-(+)-glucosamine
- 6 N-Acetyl-D-(+)-galactosamine

