


Benzodiazepines on ZB-1MS**Column:** Zebtron[™] ZB-1MS, GC Cap. Column 12 m x 0.20 mm x 0.33 μ m, Ea**Phase:** 100% Dimethylpolysiloxane**Dimensions:** 12 meters x 0.2 mm x 0.33 μ m**Order No:** 7DE-G011-14**Oven Profile:** 150°C to 200°C at 30°C/min to 300 at 15°C/min for 2 min**Carrier Gas:** Constant Flow Helium, 1.2 mL/min**Injection:** Splitless :1 3 μ L @ 280°C**Detection:** Mass Selective (MSD) (0°C)

15782

**Products used in this application:****ANALYTES:**

- 1 Oxazepam
- 2 Lorazepam
- 3 Diazepam
- 4 Nordazepam
- 5 Temazepam
- 6 Flurazepam

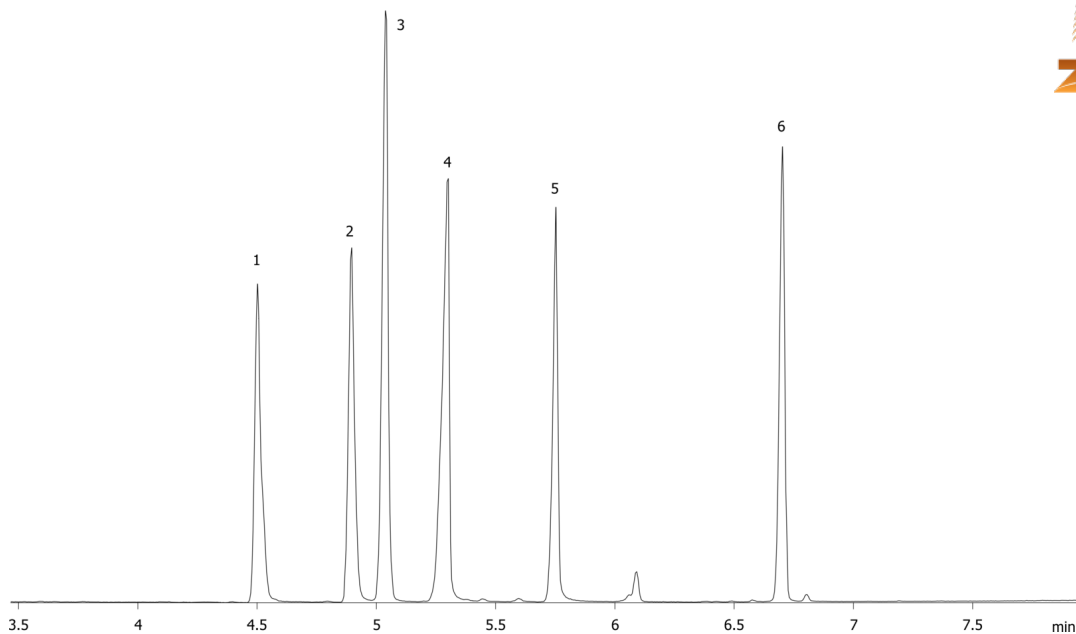