

EPA method 8081: Organochlorine PesticidesI on ZB-1

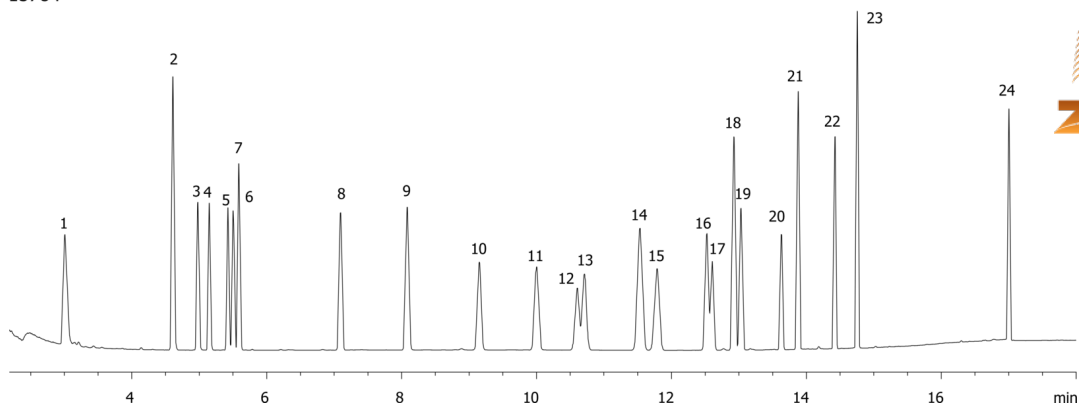
Column: Zebtron™ ZB-1, GC Cap. Column 30 m x 0.32 mm x 0.50 µm, Ea
Phase: 100% Dimethylpolysiloxane
Dimensions: 30 meters x 0.32 mm x 0.5 µm
Order No: 7HM-G001-17
Oven Profile: 140°C to 215°C at 20°C/min for 8 min to 320°C at 25°C/min for 4 min

Carrier Gas: Constant Flow Helium, 12.5 mL/min

Injection: Splitless :1 1 µL @ 250°C

Detection: Electron Capture (ECD) (330°C)

15784

**Products used in this application:**

EPA method 8081: Organochlorine PesticidesI on ZB-1

ANALYTES:

- 1 1-Bromo-2-nitrobenzene
- 2 Tetrachloro-m-xylene
- 3 alpha-BHC
- 4 Pentachloronitrobenzene
- 5 gama-BHC
- 6 Heptachlor
- 7 b-BHC
- 8 delta-BHC
- 9 Aldrin
- 10 Heptachlor epoxide
- 11 gamma-Chlordane
- 12 alpha-Chlordane
- 13 Endosulfan I
- 14 4,4'-DDE
- 15 Dieldrin
- 16 Endrin
- 17 4,4'-DDD
- 18 Endosulfan II
- 19 Endosulfan sulfate
- 20 4,4'-DDT
- 21 Endrin aldehyde
- 22 Methoxychlor
- 23 Endrin ketone
- 24 Decachlorobiphenyl

