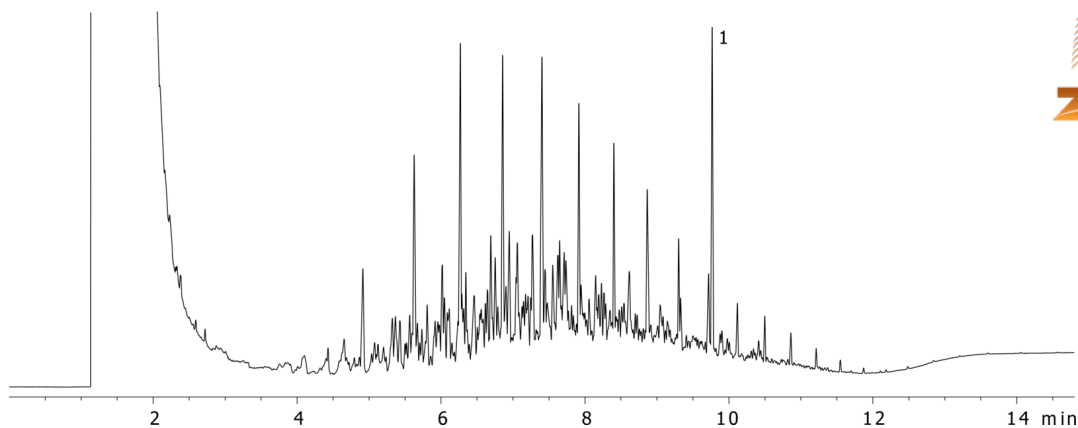


Diesel Fuel Oil #2 on ZB-5MS**Column:** ZebronTM ZB-5MS, GC Cap. Column 20 m x 0.18 mm x 0.18 μ m, Ea**Phase:** 5% Phenyl-Arylene 95% Dimethylpolysiloxane**Dimensions:** 20 meters x 0.18 mm x 0.18 μ m**Order No:** 7FD-G010-08**Oven Profile:** 50°C for 2 min to 320°C at 25°C/min for 2 min**Carrier Gas:** Constant Flow Helium, 0.8 mL/min**Injection:** Splitless :1 1 μ L @ 320°C**Detection:** Flame Ionization (FID) (335°C)

15786

**Products used in this application:****ANALYTES:****1** o-Terphenyl