

Clarity Cartridge Batch Reproducibility

Column: Clarity® 5 µm Oligo-RP, LC Column 50 x 4.6 mm, Ea

Dimensions: 50 x 4.6 mm ID

Order No: 00B-4442-E0

Elution Type: Gradient

Eluent A: 50mM TEAA, pH 7.5 / 5% ACN

Eluent B: Methanol

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0.5	90	10
	2	20	40	60

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 260 nm (ambient)


Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Clarity Oligo-RP 4 x 3.0mm, 10/Pk Part No.: AJ0-8135

- Holder Part No.: KJ0-4282



Products used in this application:

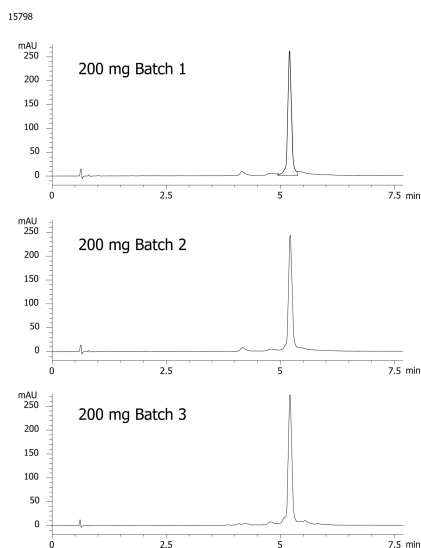

ANALYTES:

1 Purified DNA Primer

Retention Time: 6 min

