


p-Xylene Purity on ZB-WAX Plus 30x0.25x0.25**Column:** Zebron™ ZB-WAXplus™, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** Polyethylene Glycol (PEG)**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G013-11**Oven Profile:** 70°C for 6 min to 210°C at 25°C/min for 4 min**Carrier Gas:** Constant Flow Helium, 1.2 mL/min**Injection:** Split 75:1 0.1 µL @ 220°C**Detection:** Flame Ionization (FID) (250°C)

15810

**Products used in this application:****ANALYTES:**

- 1 Toluene
- 2 Ethylbenzene
- 3 p-Xylene
- 4 m-Xylene
- 5 o-Xylene
- 6 Methylbenzaldehyde

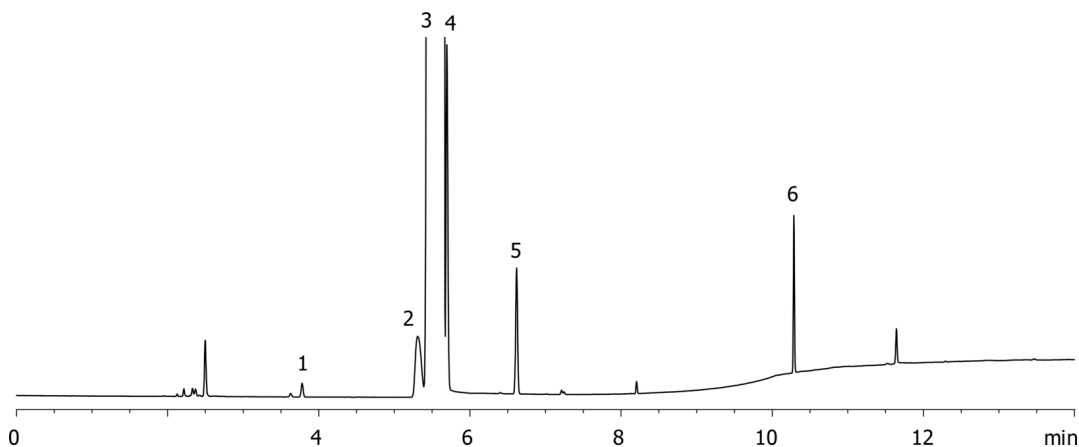