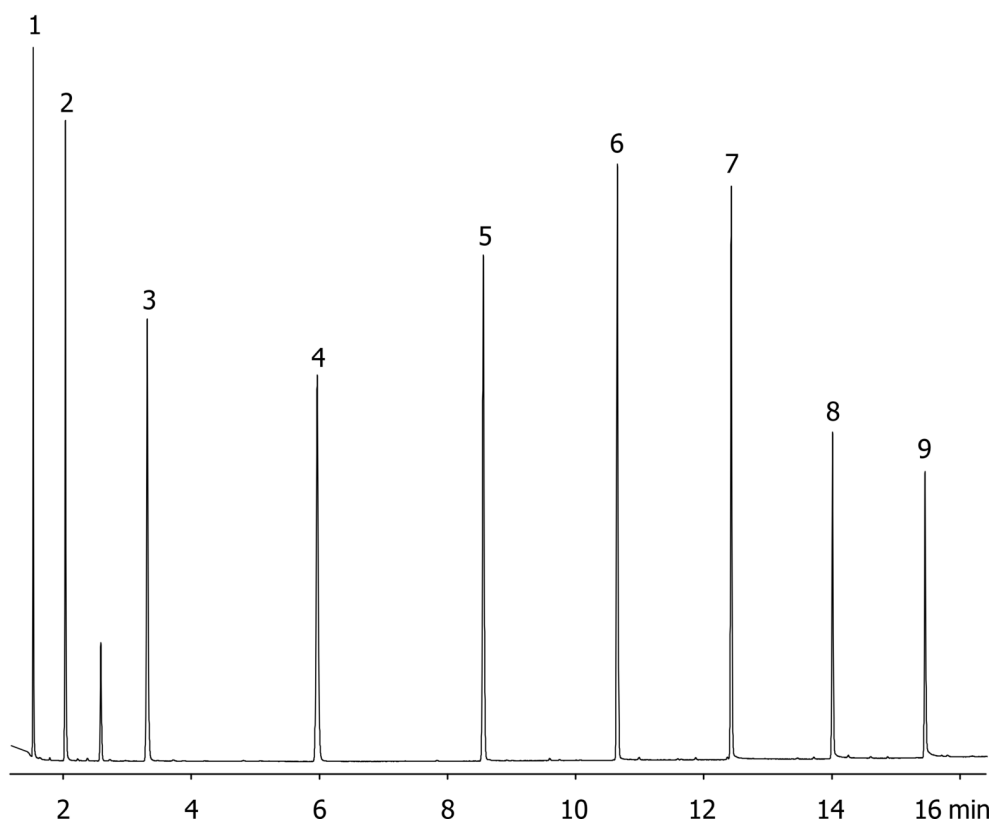


Aldehydes on ZB-WAX Plus**Column:** Zebron[™] ZB-WAXplus[™], GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** Polyethylene Glycol (PEG)**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G013-11**Oven Profile:** 40°C for 5 min to 200°C at 10°C/min**Carrier Gas:** Constant Flow Hydrogen, 1 mL/min**Injection:** Split 100:1 1 μ L @ 250°C**Detection:** Flame Ionization (FID) (225°C)

15814

**Products used in this application:****ANALYTES:**

- 1 Propanal
- 2 Butanal
- 3 Pentanal
- 4 Hexanal
- 5 Heptanal
- 6 Octanal
- 7 Nonanal
- 8 Decanal
- 9 Undecanal

©2005 Phenomenex. All rights reserved.

For more information contact your Phenomenex Representative at info@phenomenex.com

Phenomenex products are available worldwide.

www.phenomenex.com.cninfo@phenomenex.com