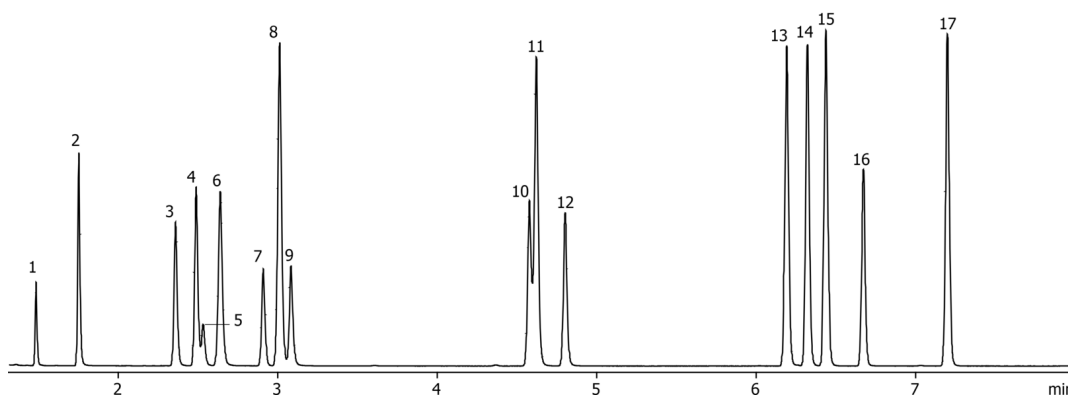


ASTM D5135-95: Styrene Impurities on ZB-WAXplus**Column:** Zebron[™] ZB-WAXplus[™], GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** Polyethylene Glycol (PEG)**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G013-11**Oven Profile:** 35°C for 2.5 min to 85°C at 10°C/min for 1 min**Carrier Gas:** Constant Flow Hydrogen, 1 mL/min**Injection:** Split 100:1 1 μ L @ 250°C**Detection:** Flame Ionization (FID) (225°C)

15815

**Products used in this application:**

ASTM D5135-95: Styrene Impurities on ZB-WAXplus

ANALYTES:

- 1 Methyl formate
- 2 Acetone
- 3 Ethyl acetate
- 4 Methyl ethyl ketone
- 5 Methanol
- 6 2-Methyl-2-propanol
- 7 Methylene chloride
- 8 Benzene
- 9 Ethanol
- 10 2-Butanol
- 11 Toluene
- 12 n-Propanol
- 13 Ethyl benzene
- 14 p-Xylene
- 15 m-Xylene
- 16 1-Butanol
- 17 o-Xylene

