

Antihistamines at Neutral pH on Gemini 5um C6-Phenyl

Column: Gemini® 5 µm C6-Phenyl 110 Å, LC Column 150 x 4.6 mm, Ea**Dimensions:** 150 x 4.6 mm ID**Order No:** 00F-4444-E0**Elution Type:** Isocratic**Eluent A:** 20mM Phosphate buffer pH7.0**Eluent B:** Methanol

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	70	30


Flow Rate: 1.5 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 210 nm (ambient)**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Gemini C6-Phenyl 4 x 3.0mm, 10/Pk Part No.: AJ0-7915

- Holder Part No.: KJ0-4282



Products used in this application:

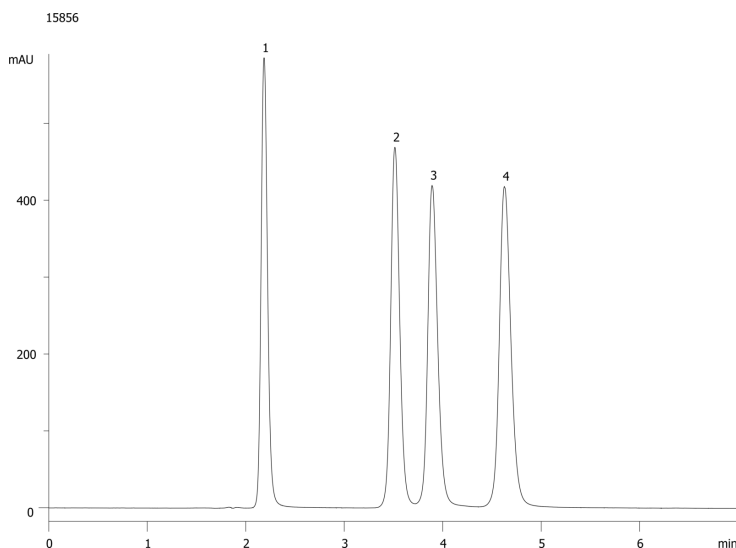

ANALYTES:

- 1 Phenylephrine
- 2 Phenylpropanolamine
- 3 Pseudoephedrine
- 4 Methylephedrine

