

European Red List Organophosphate Pesticides on Zebron ZB-1701P 30m x 0.25mm x 0.25um

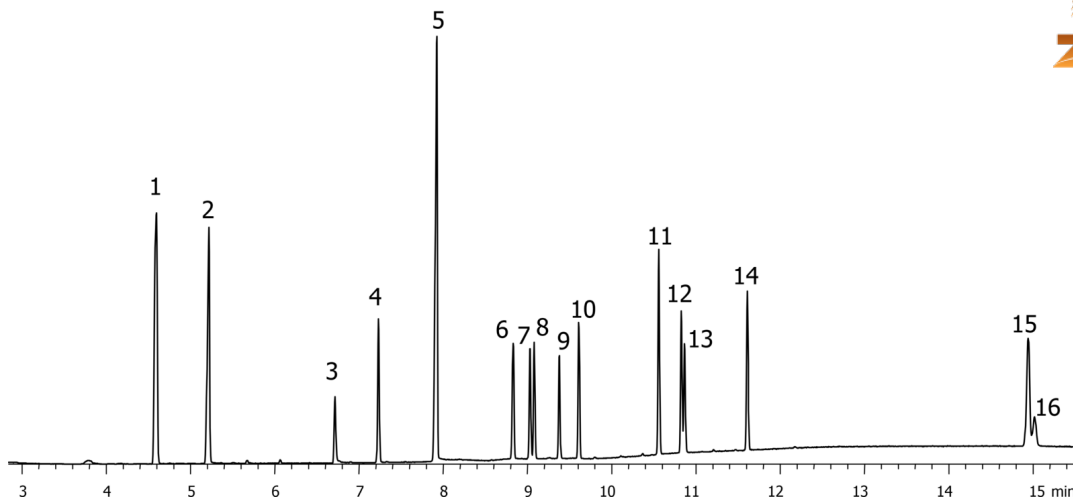
Column: Zebron[™] ZB-1701P, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea
Phase: 14% Cyanopropylphenyl 86% Dimethylpolysiloxane
Dimensions: 30 meters x 0.25 mm x 0.25 μ m
Order No: 7HG-G012-11
Oven Profile: 100°C to 220°C at 20°C/min for 1 min to 280°C at 15°C/min for 4.5 min

Carrier Gas: Constant Flow Helium, 0.8 mL/min

Injection: Splitless :1 1 μ L @ 220°C

Detection: Mass Selective (MSD) (0°C)

15864

**Products used in this application:**

European Red List Organophosphate Pesticides on Zebron ZB-1701P 30m x 0.25mm x 0.25um

ANALYTES:

- 1 Dichlorvos
- 2 Methamidophos
- 3 Acephate
- 4 Demeton-S-methyl
- 5 Omethoate
- 6 Dimethoate
- 7 Tolclofos-methyl
- 8 Pirimiphos methyl
- 9 Chlorpyrifos (Dursban)
- 10 Malathion
- 11 Prothiophos
- 12 Profenofos
- 13 Methidathion
- 14 Ethion
- 15 Pyrazophos
- 16 Azinphos-methyl

