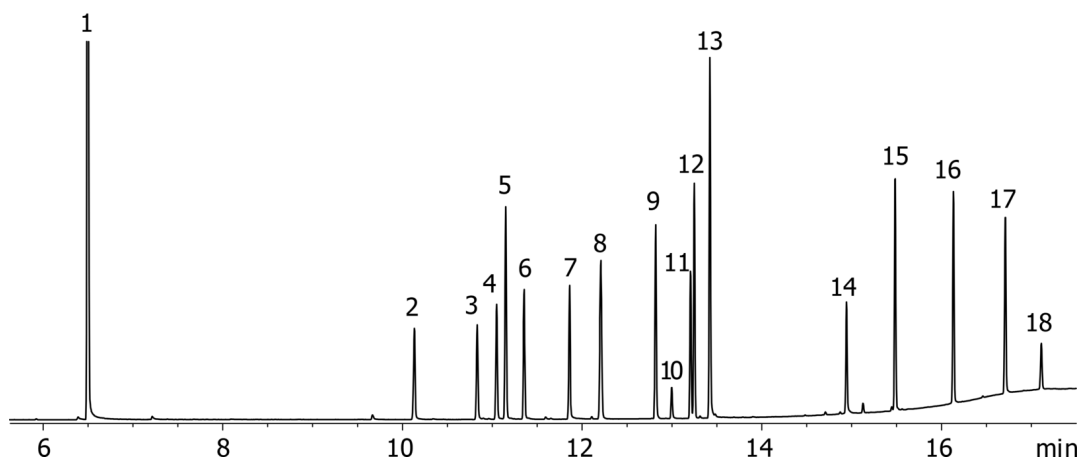


Drug Screen on Zebron ZB-5MS 20m x 0.18mm x 0.18um

Column: Zebron[™] ZB-5MS, GC Cap. Column 20 m x 0.18 mm x 0.18 μ m, Ea
Phase: 5% Phenyl-Arylene 95% Dimethylpolysiloxane
Dimensions: 20 meters x 0.18 mm x 0.18 μ m
Order No: 7FD-G010-08
Oven Profile: 50°C for 1 min to 125°C at 25°C/min to 200°C at 12°C/min to 325°C at 20°C/min for 1 min
Carrier Gas: Constant Flow Helium, 1 mL/min
Injection: Split 20:1 1 μ L @ 230°C
Detection: Flame Ionization (FID) (330°C)
15868

**Products used in this application:**

Drug Screen on Zebron ZB-5MS 20m x 0.18mm x 0.18um

ANALYTES:

- 1** Nicotine
- 2** Pentobarbital
- 3** Caffeine
- 4** Glutethimide
- 5** Lidocaine
- 6** Phencyclidine
- 7** Phenobarbital
- 8** EDDP
- 9** Methadone
- 10** Methaqualone
- 11** Cocaine
- 12** Trimipramine
- 13** Desipramine
- 14** Heroin
- 15** Fentanyl
- 16** Ibogaine
- 17** Triazolam
- 18** LSD

