


Solvent Mix on Zebron ZB-1MS 30m x 0.25mm x 1.0um**Column:** Zebron[™] ZB-1MS, GC Cap. Column 30 m x 0.25 mm x 1.00 μ m, Ea**Phase:** 100% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 1 μ m**Order No:** 7HG-G011-22**Oven Profile:** 50°C for 4.5min to 120°C at 20°C/min**Carrier Gas:** Constant Flow Helium, 3.2 mL/min**Injection:** Headspace 15:1 200 μ L @ 120°C**Detection:** Flame Ionization (FID) (250°C)**Analyst Note:** 200 μ L Headspace injection at 50°C Split 15:1.
Analytes are 0.05mg/20mL.**Products used in this application:**

Solvent Mix on Zebron ZB-1MS 30m x 0.25mm x 1.0um

ANALYTES:

- 1 Methanol
- 2 Ethanol
- 3 Acetone
- 4 Isopropanol
- 5 n-Propanol
- 6 2-Butanone
- 7 Ethyl acetate
- 8 Isopropyl acetate
- 9 Methoxypropanol
- 10 n-Propyl acetate
- 11 4-Methyl-2-pentanone
- 12 Ethoxypropanol
- 13 Toluene
- 14 2,4-Pentanedione
- 15 n-Butyl acetate

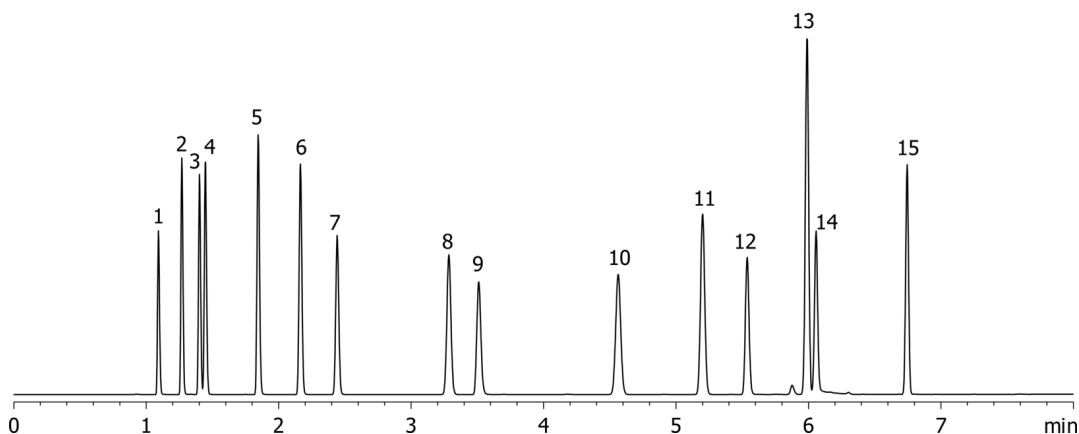