


Blood Alcohols on ZB-WAXplus**Column:** ZebronTM ZB-WAXplusTM, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** Polyethylene Glycol (PEG)**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G013-11**Oven Profile:** 65°C isothermal**Carrier Gas:** Constant Flow Helium, 1.5 mL/min**Injection:** Split 20:1 50 μ L @ 225°C**Detection:** Flame Ionization (FID) (250°C)

15880

**Products used in this application:****ANALYTES:**

- 1 Acetaldehyde
- 2 Acetone
- 3 Methanol
- 4 Isopropanol
- 5 Ethanol
- 6 n-Propanol

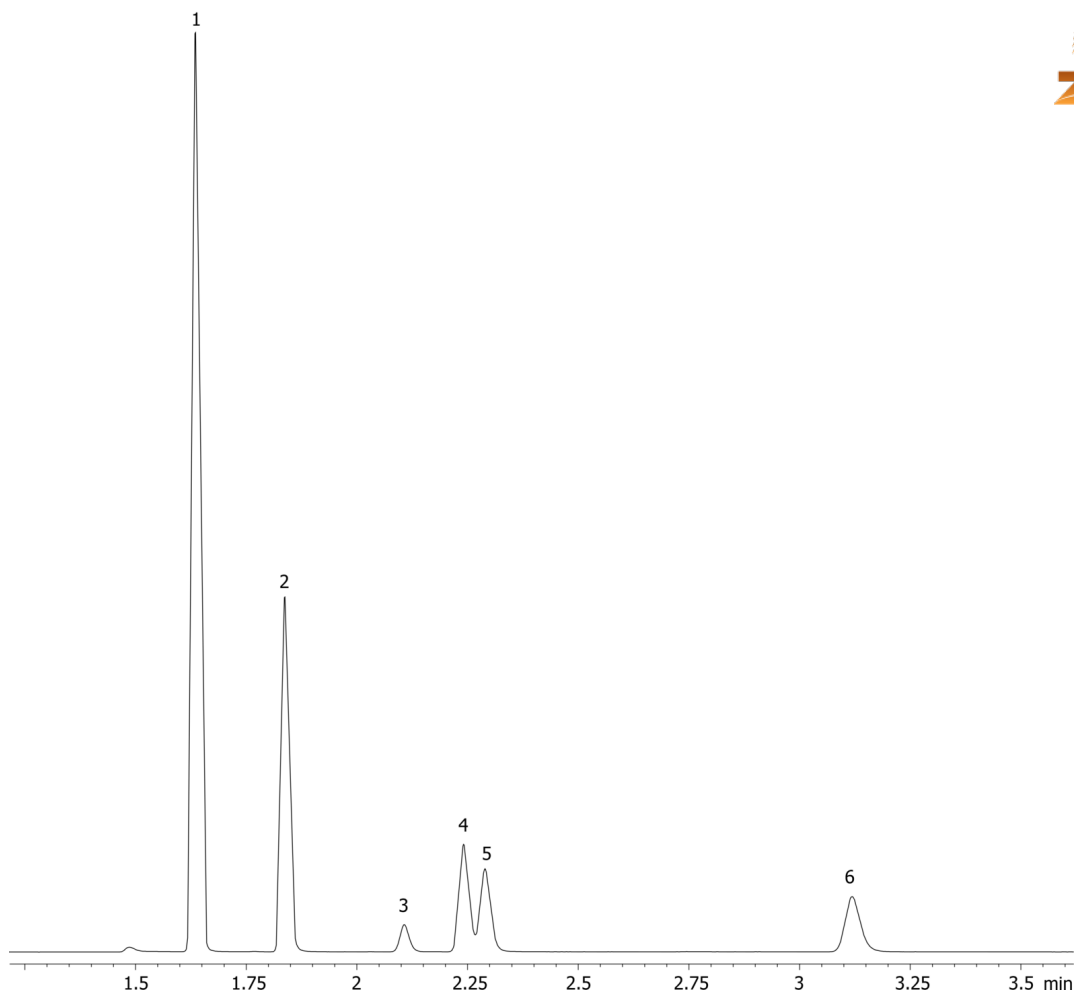