

Soft Drink Application using Gemini 5u C6-Phenyl

Column: Gemini® 5 µm C6-Phenyl 110 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4444-E0

Elution Type: Isocratic

Eluent A: 20mM Potassium Phosphate pH 2.5

Eluent B: Methanol

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	55	45

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Gemini C6-Phenyl 4 x 3.0mm, 10/Pk Part No.: AJ0-7915


- Holder Part No.: KJ0-4282



Products used in this application:



15891

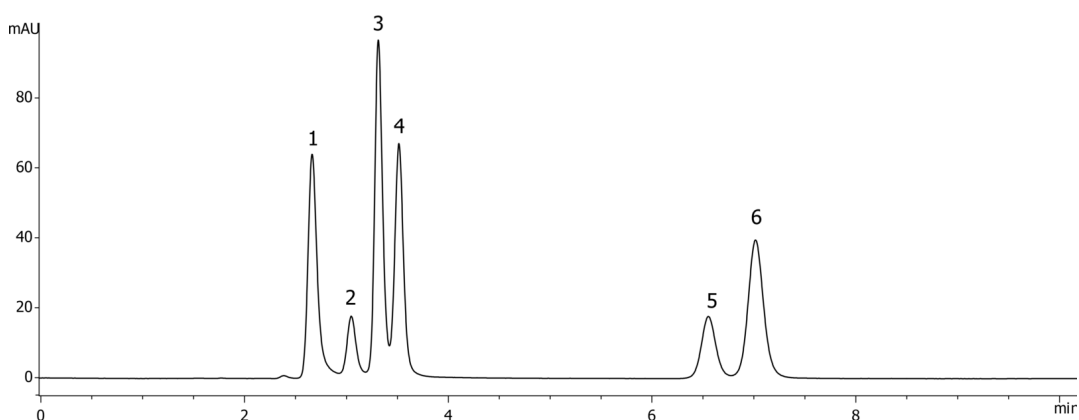

ANALYTES:

- 1 Quinine
- 2 Aspartame
- 3 Saccharin
- 4 Caffeine
- 5 Sorbic acid
- 6 Benzoic acid

