

Sulfa Drugs using Luna 5u Phenyl-Hexyl

Column: Luna® 5 µm Phenyl-Hexyl 100 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4257-E0

Elution Type: Isocratic

Eluent A: 0.1% Formic Acid in Water

Eluent B: Methanol

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
1	1	0	70	30

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Phenyl 4 x 3.0mm, 10/Pk Part No.: AJ0-4351


- Holder Part No.: KJ0-4282



Products used in this application:



15905

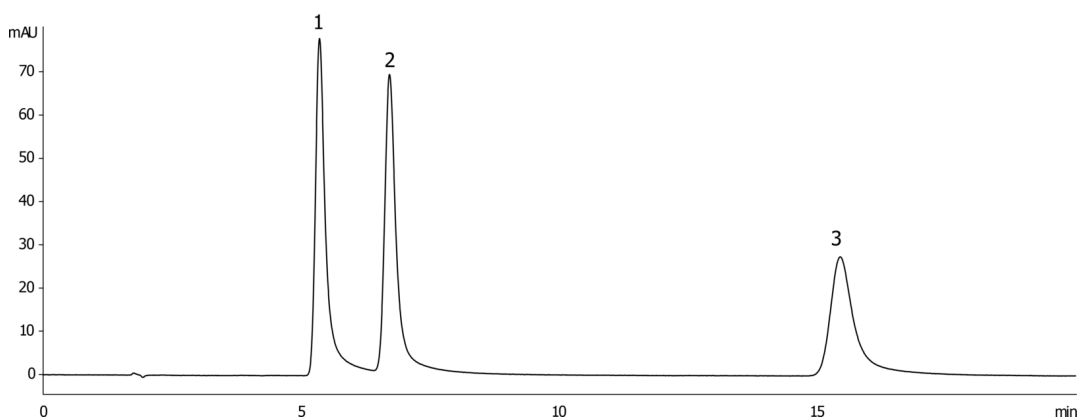

ANALYTES:

- 1 Sulfathiazole
- 2 Sulfamerazine
- 3 Sulfamethoxazole

