


Common Solvents on ZB-WAXplus**Column:** Zebtron™ ZB-WAXplus™, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** Polyethylene Glycol (PEG)**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G013-11**Oven Profile:** 50°C for 4 min to 140°C at 20°C/min**Carrier Gas:** Constant Flow Helium, 1.4 mL/min**Injection:** Split 25:1 4 µL @ 100°C**Detection:** Mass Selective (MSD) (250°C)**Analyst Note:** 4 µL Headspace split 25:1 at 100°C**Products used in this application:**

Common Solvents on ZB-WAXplus

ANALYTES:

- 1 Acetone
- 2 Ethyl acetate
- 3 Isopropyl acetate
- 4 Methanol
- 5 Methyl ethyl ketone
- 6 Isopropyl alcohol
- 7 Ethanol
- 8 n-Propyl acetate
- 9 Methyl isobutyl ketone
- 10 n-Propanol
- 11 Toluene
- 12 n-Butyl acetate
- 13 Methoxypropanol (Propylene glycol monomethyl ether)
- 14 3-Ethoxy propanol (Propylene glycol Ä-monoethyl ether)

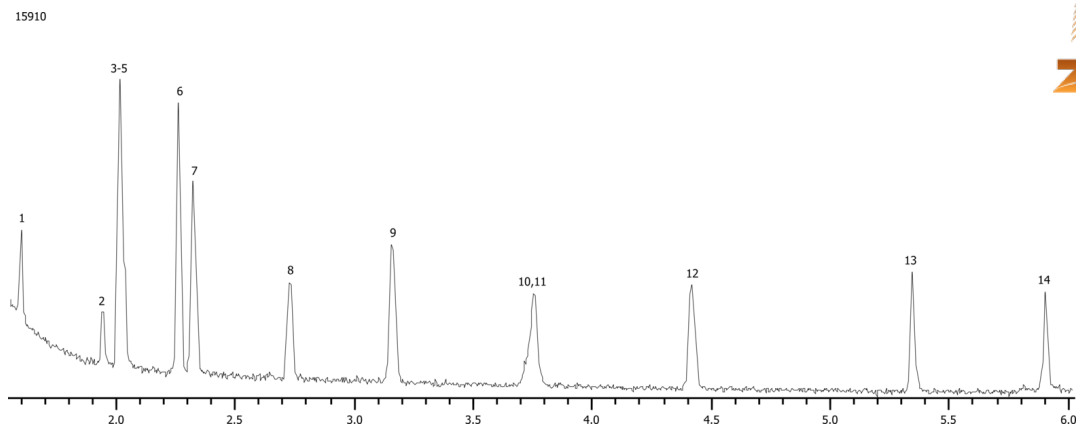