


Glycols on ZB-WAXplus**Column:** Zebron[™] ZB-WAXplus[™], GC Cap. Column 30 m x 0.32 mm x 0.50 μ m, Ea**Phase:** Polyethylene Glycol (PEG)**Dimensions:** 30 meters x 0.32 mm x 0.5 μ m**Order No:** 7HM-G013-17**Oven Profile:** 80°C to 210°C at 25°C/min for 6 min**Carrier Gas:** Constant Flow Helium, 1.8 mL/min**Injection:** Split 20:1 1 μ L @ 220°C**Detection:** Flame Ionization (FID) (250°C)**Analyst Note:** Analytes are 100ppm in Dichloromethane**Products used in this application:****ANALYTES:**

- 1 Ethylene glycol dimethyl ether
- 2 Propylene glycol
- 3 Ethylene glycol
- 4 Ethylene glycol dimethacrylate
- 5 Ethylene glycol diglycidyl ether
- 6 Ethylene glycol phenyl methacrylate

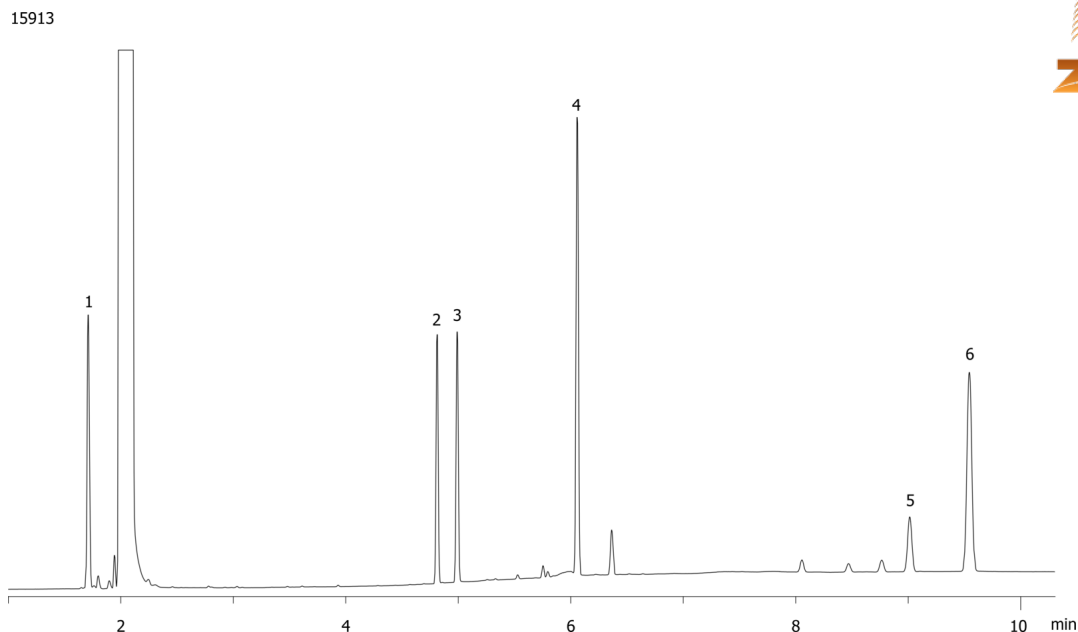