

Poly dT DNA Oligo-RP Purification

Column: Clarity® 3 µm Oligo-RP, LC Column 50 x 4.6 mm, Ea

Dimensions: 50 x 4.6 mm ID

Order No: 00B-4441-E0

Elution Type: Gradient

Eluent A: 50mM TEAA/5% Acetonitrile

Eluent B: Methanol

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	90	10
	2	5	85	15
	3	35	70	30

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 260 nm (ambient)


Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Clarity Oligo-RP 4 x 3.0mm, 10/Pk Part No.: AJ0-8135

- Holder Part No.: KJ0-4282



Products used in this application:

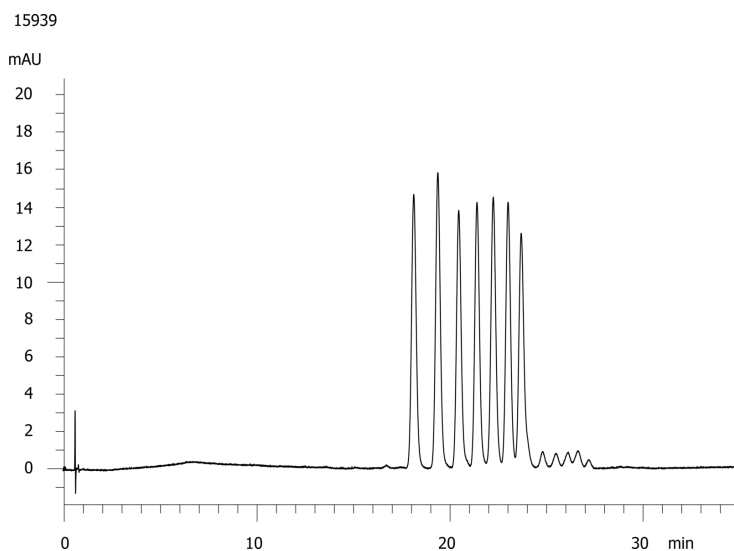

ANALYTES:

- 1 Poly-dT 12
- 2 Poly-dT 13
- 3 Poly-dT 14
- 4 Poly-dT 15
- 5 Poly-dT 16
- 6 Poly-dT 17
- 7 Poly-dT 18

