

Prepartive (200µg) 20nt DNA Oligo-RP 3µ Purification

Column: Clarity® 3 µm Oligo-RP, LC Column 50 x 10 mm, Ea

Dimensions: 50 x 10 mm ID

Order No: 00B-4441-N0

Elution Type: Gradient

Eluent A: 50mM TEAA/5% Acetonitrile

Eluent B: Methanol

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	90	10
	2	20	40	60

Flow Rate: 5 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 260 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard SemiPrep Cartridges, Clarity Oligo-RP 10 x 10mm, 3/Pk Part No.: AJ0-8136


- Holder Part No.: AJ0-7220



Products used in this application:



15947

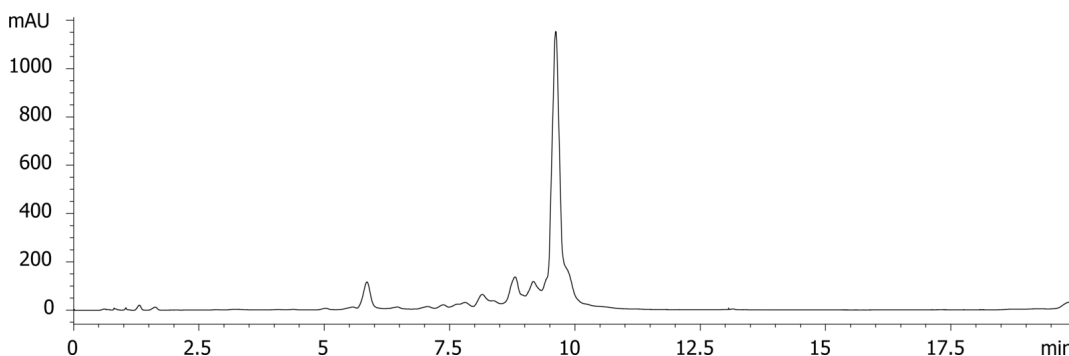

ANALYTES:

1 20nt DNA

