

## 60nt DNA Cartridge Purification - Crude, Load, Wash, Elution Fractions

**Column:** Clarity® 3 µm Oligo-RP, LC Column 50 x 4.6 mm, Ea

**Dimensions:** 50 x 4.6 mm ID

**Order No:** 00B-4441-E0

**Elution Type:** Gradient

**Eluent A:** 50mM TEAA/5% Acetonitrile

**Eluent B:** Methanol

| Gradient Profile: | Step No. | Time (min) | Pct A | Pct B |
|-------------------|----------|------------|-------|-------|
|                   | 1        | 0          | 90    | 10    |
|                   | 2        | 20         | 40    | 60    |

**Flow Rate:** 1 mL/min

**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 260 nm (ambient)


**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Clarity Oligo-RP 4 x 3.0mm, 10/Pk Part No.: AJ0-8135

- Holder Part No.: KJ0-4282



Products used in this application:

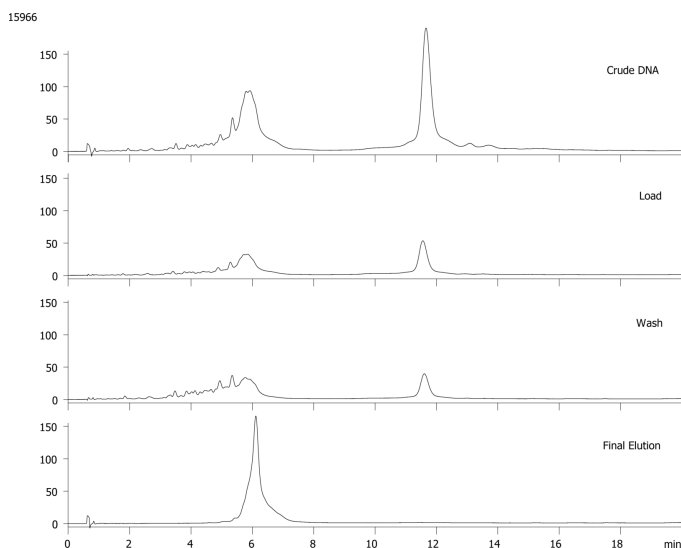

### ANALYTES:

1 DNA (60nt)

