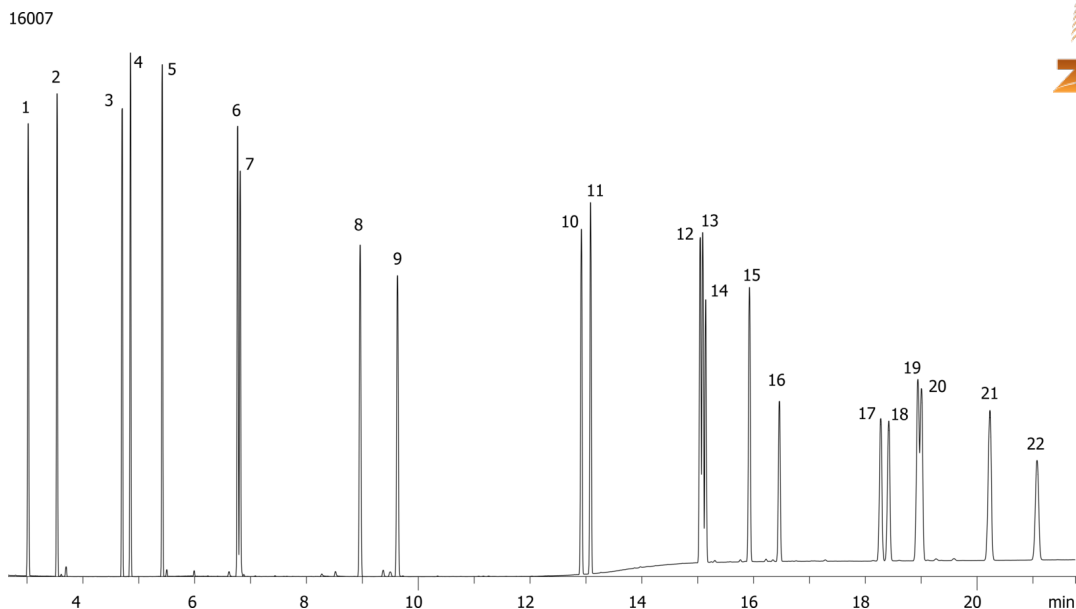


Extended PAHs on ZB-50 60x0.25x0.25

Column: Zebron[™] ZB-50, GC Cap. Column 60 m x 0.25 mm x 0.25 μ m, Ea
Phase: 50% Phenyl 50% Dimethylpolysiloxane
Dimensions: 60 meters x 0.25 mm x 0.25 μ m
Order No: 7KG-G004-11
Oven Profile: 140°C to 240°C at 20°C/min to 295°C at 10°C/min to 320°C at 30°C/min for 43 min
Carrier Gas: Constant Flow Helium, 1.8 mL/min
Injection: Split 20:1 1 μ L @ 320°C
Detection: Flame Ionization (FID) (325°C)
Analyst Note: All samples are 20ppm in Dichloromethane

**Products used in this application:**

Extended PAHs on ZB-50 60x0.25x0.25

ANALYTES:

- 1 Naphthalene
- 2 2-Methylnaphthalene
- 3 Acenaphthalene
- 4 Acenaphthene
- 5 Fluorene
- 6 Phenanthrene
- 7 Anthracene
- 8 Fluoranthene
- 9 Pyrene
- 10 Benz[a]anthracene
- 11 Chrysene
- 12 Benzo[b]fluoranthene
- 13 Benzo[k]fluoranthene
- 14 Benzo[j]fluoranthene
- 15 Benzo[a]pyrene
- 16 3-Methylcholanthrene
- 17 Dibenzo[a,h]acridine
- 18 Dibenzo[a,,j]acridine
- 19 Dibenzo[a,h]anthracene
- 20 Indeno[1,2,3-cd]pyrene
- 21 Benzo[g,h,i]perylene
- 22 7H-Dibenzo[c,g]carbazole

