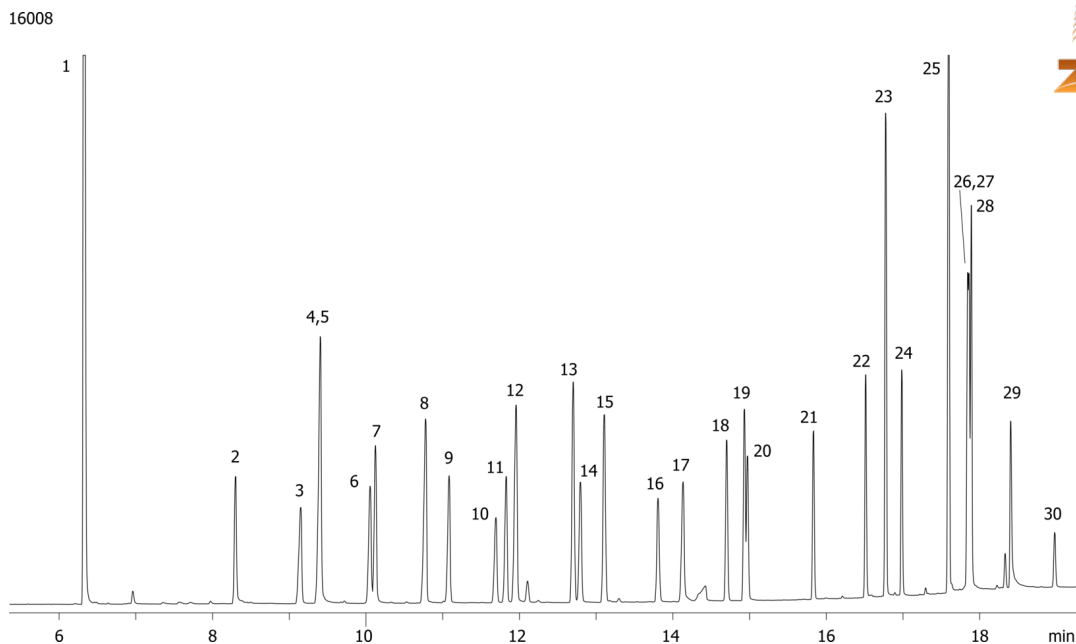


Common Drug Screen on ZB-1MS

Column: Zebtron™ ZB-1MS, GC Cap. Column 20 m x 0.18 mm x 0.18 µm, Ea
Phase: 100% Dimethylpolysiloxane
Dimensions: 20 meters x 0.18 mm x 0.18 µm
Order No: 7FD-G011-08
Oven Profile: 80°C for 1 min to 150°C at 15°C/min to 200°C at 6°C/min to 275°C at 20°C/min for 2 min
Carrier Gas: Constant Flow Helium, 0.9 mL/min
Injection: Split 10:1 1.5 µL @ 225°C
Detection: Flame Ionization (FID) (250°C)
Analyst Note: Analytes are 25ppm except nicotine at 100ppm

**Products used in this application:**

Common Drug Screen on ZB-1MS

ANALYTES:

- 1 Nicotine
- 2 Acetophenone
- 3 Allobarbitol
- 4 Aprobarbitol
- 5 Ibuprofen
- 6 Butalbital
- 7 Phenacetin
- 8 Amobarbital
- 9 Pentobarbital
- 10 Meprobamate
- 11 Secobarbital
- 12 Caffeine
- 13 Benzphetamine
- 14 Hexobarbital
- 15 Dimenhydrinate
- 16 Doxylamine
- 17 Phenobarbital
- 18 Methapyrilene
- 19 Chlorpheniramine
- 20 Procaine
- 21 Brompheniramine
- 22 Cocaine
- 23 Chlorcyclizine
- 24 Benactyzine
- 25 Codeine
- 26 Morphine
- 27 Diazepam
- 28 Hydrocodone
- 29 Oxymorphone
- 30 Heroin

