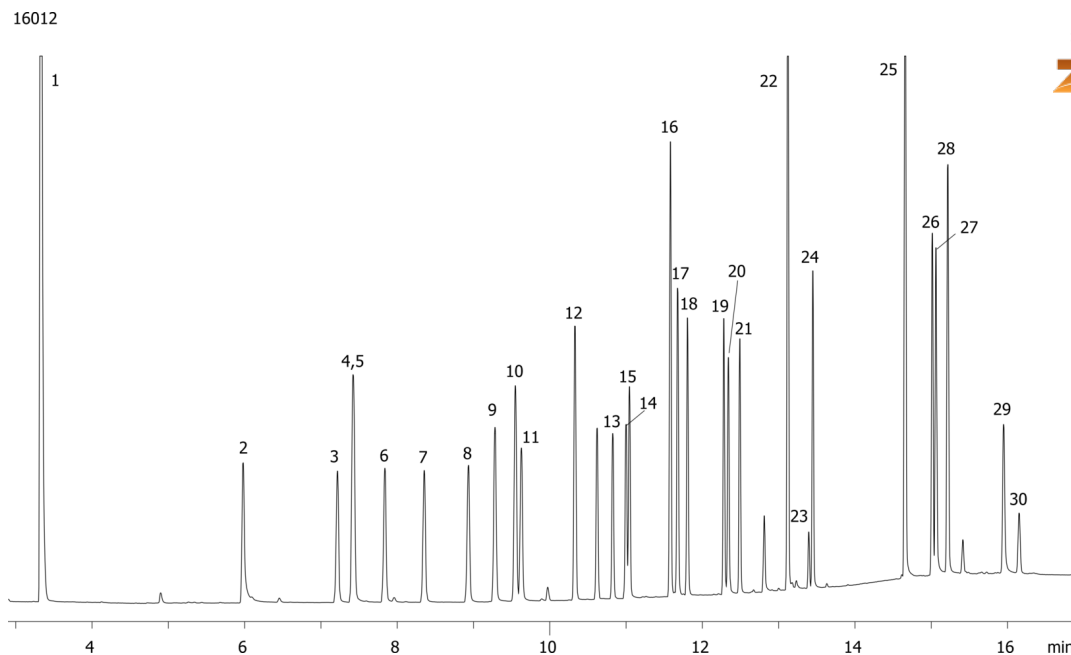


Common Drug Screen on ZB-50

Column: Zebron[™] ZB-50, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea
Phase: 50% Phenyl 50% Dimethylpolysiloxane
Dimensions: 30 meters x 0.25 mm x 0.25 μ m
Order No: 7HG-G004-11
Oven Profile: 150°C for 1min to 200°C at 6°C/min to 300°C at 20°C/min for 4min

Carrier Gas: Constant Flow Helium, 2.2 mL/min
Injection: Split 10:1 1.5 μ L @ 275°C
Detection: Flame Ionization (FID) (300°C)
Analyst Note: All analytes are 25ppm except nicotine at 100ppm

**Products used in this application:**

Common Drug Screen on ZB-50

ANALYTES:

- 1 Nicotine
- 2 Ibuprofen
- 3 Allobarbitol
- 4 Acetomenophen
- 5 Aprobital
- 6 Butalbitol
- 7 Amobarbitol
- 8 Pentobarbitol
- 9 Phenacetin
- 10 Benzphetamine
- 11 Secobarbitol
- 12 Dimenhydrinate
- 13 Meprobamate
- 14 Doxylamine
- 15 Hexobarbitol
- 16 Caffeine
- 17 Chlorpheniramine
- 18 Methapyrilene
- 19 Phenobarbitol
- 20 Procaine
- 21 Brompheniramine
- 22 Chlorcyclizine
- 23 Cocaine
- 24 Benactyzine
- 25 Codeine
- 26 Diazepam
- 27 Morphine
- 28 Hydrocodone
- 29 Oxymorphone
- 30 Heroin

