


Cold Pressed Orange Oil on ZB-1MS**Column:** Zebron[™] ZB-1MS, GC Cap. Column 10 m x 0.10 mm x 0.10 μ m, Ea**Phase:** 100% Dimethylpolysiloxane**Dimensions:** 10 meters x 0.1 mm x 0.1 μ m**Order No:** 7CB-G011-02**Oven Profile:** 60°C to 130°C at 10°C/min to 280°C at 30°C for 3min**Carrier Gas:** Constant Flow Helium, 0.3 mL/min**Injection:** Split 120:1 0.2 μ L @ 160°C**Detection:** Mass Selective (MSD) (0°C)**Analyst Note:** Sample was 10% in dichloromethane**Products used in this application:**

Cold Pressed Orange Oil on ZB-1MS

ANALYTES:

- 1 α -Pinene
- 2 Sabinene
- 3 beta-Myrcene
- 4 Octanal
- 5 3-Carene
- 6 Limonene
- 7 Nonanal
- 8 Linalool
- 9 cis-Limonene oxide
- 10 trans-Limonene oxide
- 11 Citronellal
- 12 α -Terpineol
- 13 Decanal
- 14 Carvone
- 15 Neral
- 16 Geranial
- 17 α -Cubebene
- 18 β -Cubenene
- 19 Dodecanal
- 20 Valencene
- 21 Cadinene
- 22 Nootkatone

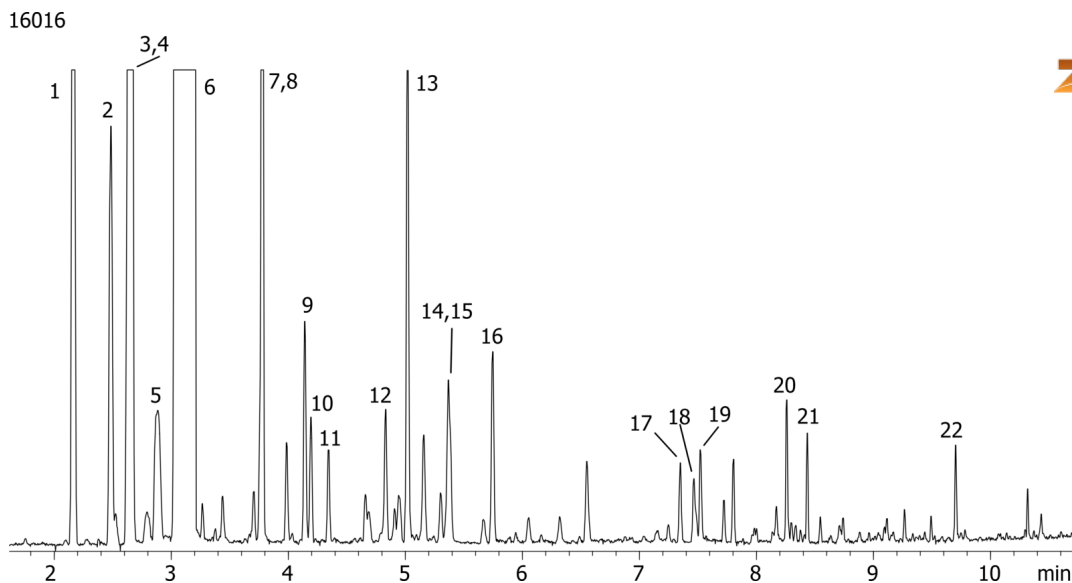