


Prodigy Test Mix using Onyx C18E 100x3.0mm (0.8mL/min)

Column: Onyx™ Monolithic C18, LC Column 100 x 3 mm, Ea
Dimensions: 100 x 3 mm ID
Order No: CH0-8158
Elution Type: Isocratic
Eluent A: Water
Eluent B: Acetonitrile
Gradient

Step No.	Time (min)	Pct A	Pct B
1	0	35	65

Profile:
Flow Rate: 0.8 mL/min
Col. Temp.: ambient
Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

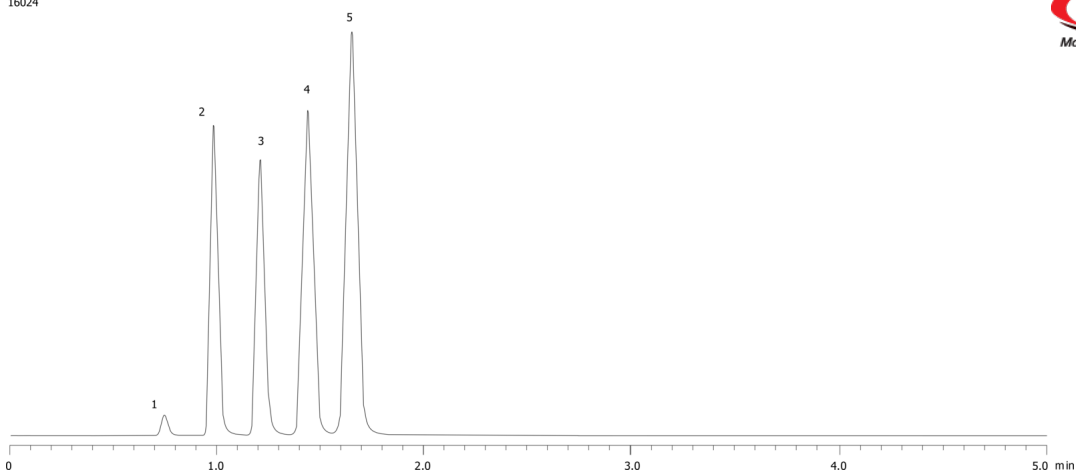

Products used in this application:



ANALYTES:

- 1 Uracil
- 2 Acetophenone
- 3 Benzene
- 4 Toluene
- 5 Naphthalene

