

MercuryMS - Synergi 2.5µm Max-RP 20 x 4.0mm for Acid, Base, Neutral

Column: Synergi™ 2.5 µm Max-RP 100 Å, LC Column 20 x 4 mm, Mercury, Ea

Dimensions: 20 x 4 mm ID

Order No: 00M-4372-D0-CE

Elution Type: Gradient

Eluent A: 0.1% Formic Acid in Water


Eluent B: 0.1% Formic Acid in Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	1	5	95
	3	1.01	95	5
	4	1.2	95	5

Flow Rate: 2.5 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

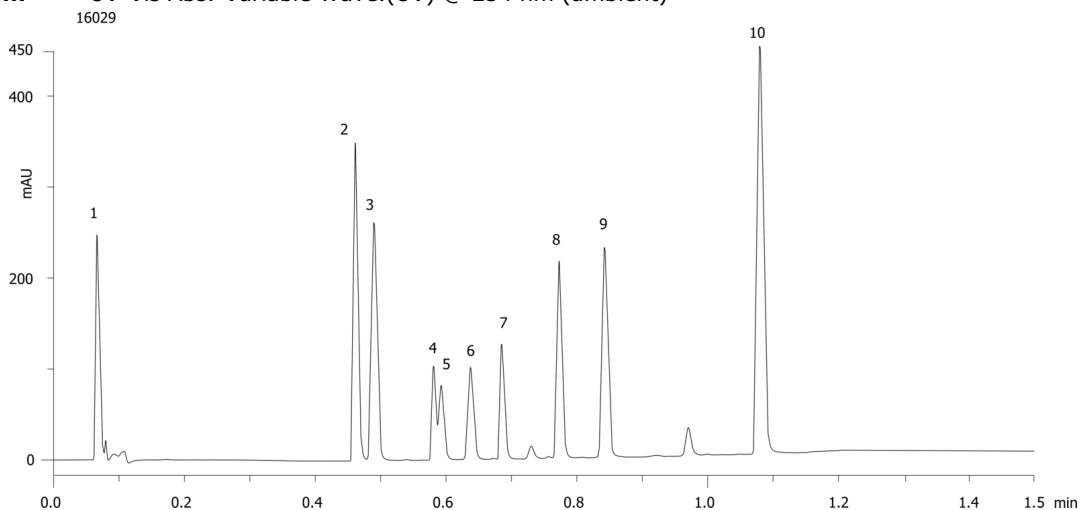

MercuryMS - Synergi 2.5 μ m Max-RP 20 x 4.0mm for Acid, Base, Neutral

ANALYTES:

- 1 Pyridine
- 2 Quinidine
- 3 Sulfathiazole
- 4 Triprolidine
- 5 Benzyl alcohol
- 6 Phenol
- 7 Nortriptyline
- 8 3-Methyl-4-nitrobenzoic acid
- 9 5-Methyl-salicylaldehyde
- 10 Hexaphenone

