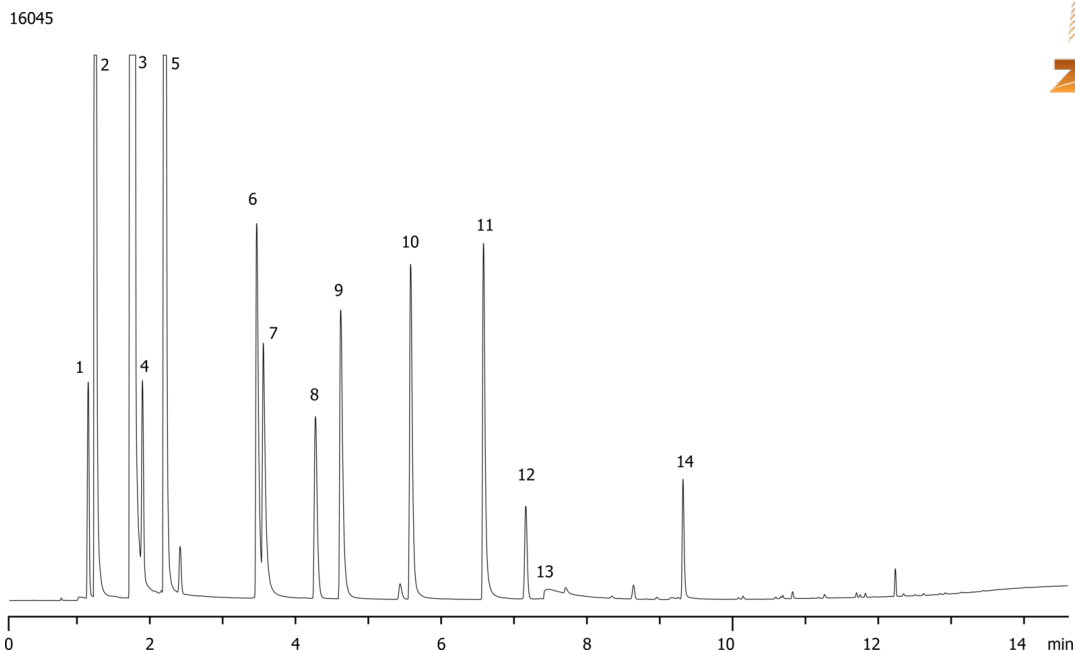


Scotch Whiskey Main Compounds on ZB-624

Column: Zebron[™] ZB-624, GC Cap. Column 20 m x 0.18 mm x 1.00 µm, Ea
Phase: 6% Cyanopropylphenyl 94% Dimethylpolysiloxane
Dimensions: 20 meters x 0.18 mm x 1 µm
Order No: 7FD-G005-22
Oven Profile: 35°C for 4 min to 85°C at 10°C/min to 200°C at 25°C/min for 1 min

Carrier Gas: Constant Flow Helium, 1.5 mL/min
Injection: Split 40:1 0.2 µL @ 140°C
Detection: Flame Ionization (FID) (200°C)
Analyst Note: Analytes are 20ppm in 60/40 water/ethanol matrix

**Products used in this application:**

Scotch Whiskey Main Compounds on ZB-624

ANALYTES:

- 1 Acetaldehyde
- 2 Methanol
- 3 Ethanol
- 4 Acrolein
- 5 Isopropanol
- 6 Allyl alcohol
- 7 1-Propanol
- 8 Ethyl acetate
- 9 sec-Butanol
- 10 Isobutanol
- 11 n-Butanol
- 12 Acetal
- 13 Acetic acid
- 14 Acetaldehyde Trimer

