


Glenfiddich® Scotch Whiskey on ZB-WAXplus**Column:** Zebron™ ZB-WAXplus™, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** Polyethylene Glycol (PEG)**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G013-11**Oven Profile:** 35°C for 5 min to 200°C at 30°C/min for 1min**Carrier Gas:** Constant Flow Helium, 1.4 mL/min**Injection:** Split 30:1 0.2 µL @ 140°C**Detection:** Flame Ionization (FID) (280°C)**Analyst Note:** Injection of pure Glenfiddich®**Products used in this application:****ANALYTES:**

- 1 Acetaldehyde
- 2 Ethyl acetate
- 3 Acetal
- 4 Methanol
- 5 Ethanol
- 6 1-Propanol
- 7 Isobutanol
- 8 Isoamyl alcohol
- 9 2-Methylbutanol

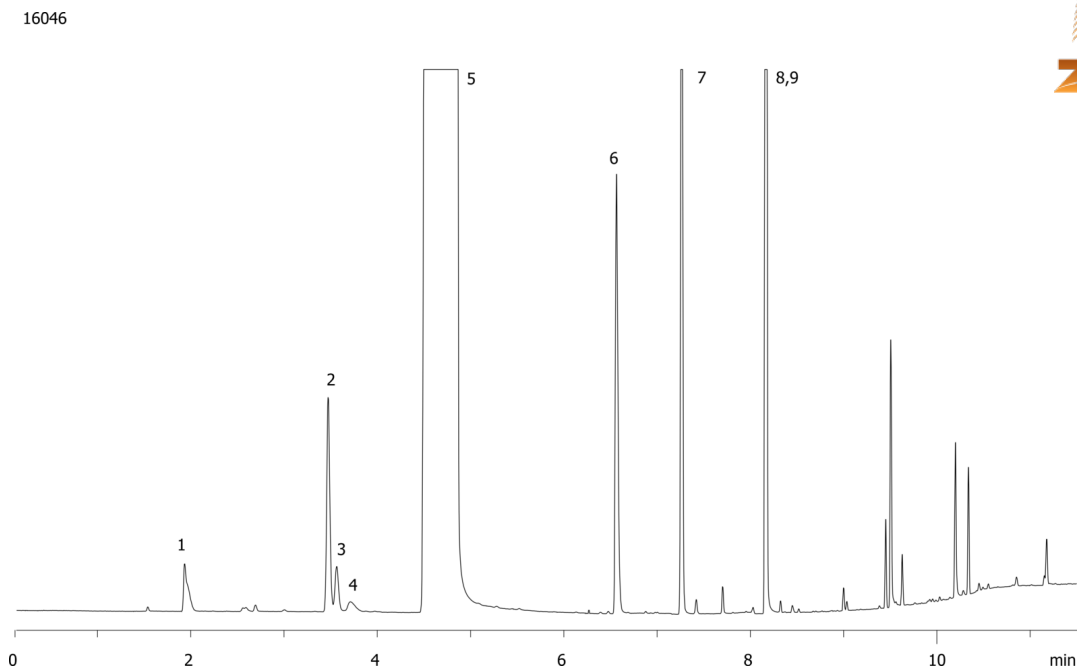