

Non-Polar Acids at pH 2.5 using Onyx C18 3.0mm ID


Column: Onyx™ Monolithic C18, LC Column 100 x 3 mm, Ea
Dimensions: 100 x 3 mm ID
Order No: CH0-8158
Elution Type: Isocratic
Eluent A: 20mM Potassium phosphate pH 2.5
Eluent B: Acetonitrile
Gradient

Step No.	Time (min)	Pct A	Pct B
1	0	50	50

Profile:
Flow Rate: 0.4 mL/min
Col. Temp.: ambient
Detection: UV-Vis Abs.-Variable Wave.(UV) @ 230 nm (ambient)



Products used in this application:

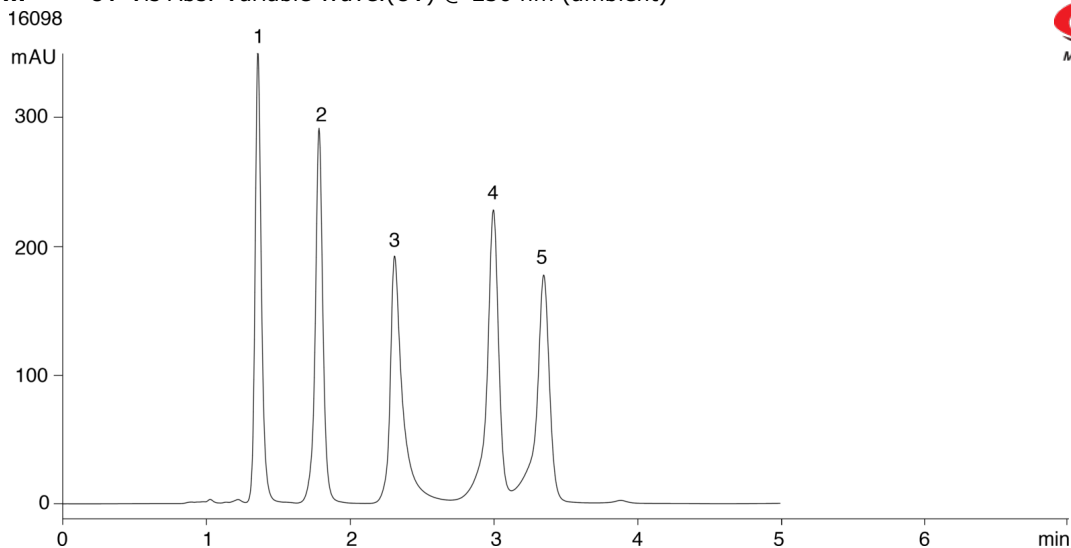

ANALYTES:

- 1 Ethyl paraben
- 2 Naproxen
- 3 Diflunisal
- 4 Indomethacin
- 5 Ibuprofen

