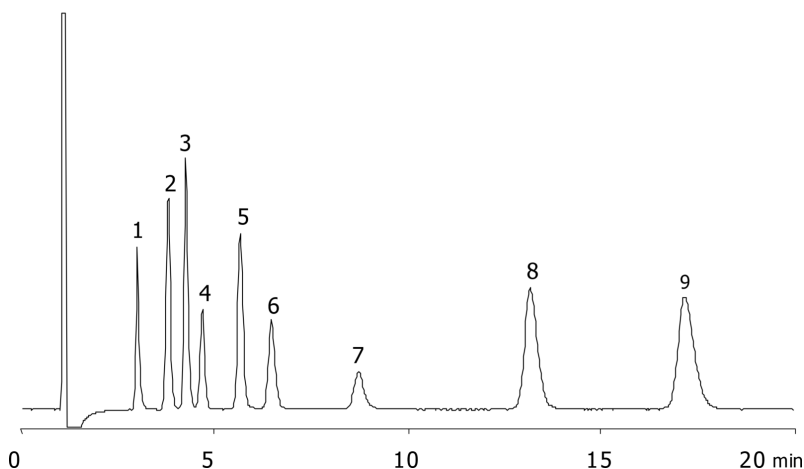


Alkyl Amines and Cations**Column:** N/A**Dimensions:** N/A**Order No:** N/A**Elution Type:** N/A**Eluent A:** 4mM Methanesulfonic acids aq.

Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100

Flow Rate: 1 mL/min**Col. Temp.:** 40 °C**Detection:** Conductivity @ (ambient)

16113

**ANALYTES:**

- 1 Li+
- 2 Na+
- 3 NH₄⁺
- 4 Methylamine
- 5 K+
- 6 Trimethylamine
- 7 Triethylamine
- 8 Mg²⁺
- 9 Ca²⁺

