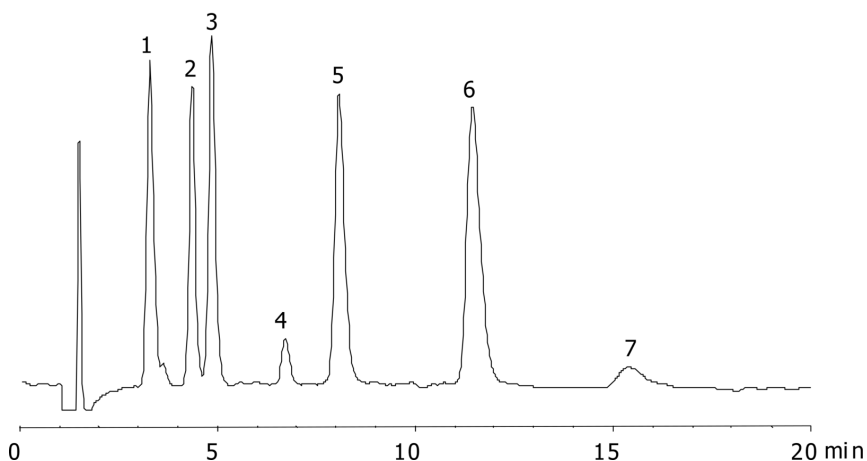


Transition Metals

Column: N/A
Dimensions: N/A
Order No: N/A
Elution Type: N/A
Eluent A: 6mM Tartaric acid + 4mM Oxalic Acid aq.
Gradient

Step No.	Time (min)	Pct A
1	0	100

Flow Rate: 1 mL/min
Col. Temp.: 40 °C
Detection: Conductivity @ (ambient)
16114

**ANALYTES:**

- 1 Ni+
- 2 Zn2+
- 3 Co2+
- 4 Fe2+
- 5 Mn2+
- 6 Cd2+
- 7 Pb2+

