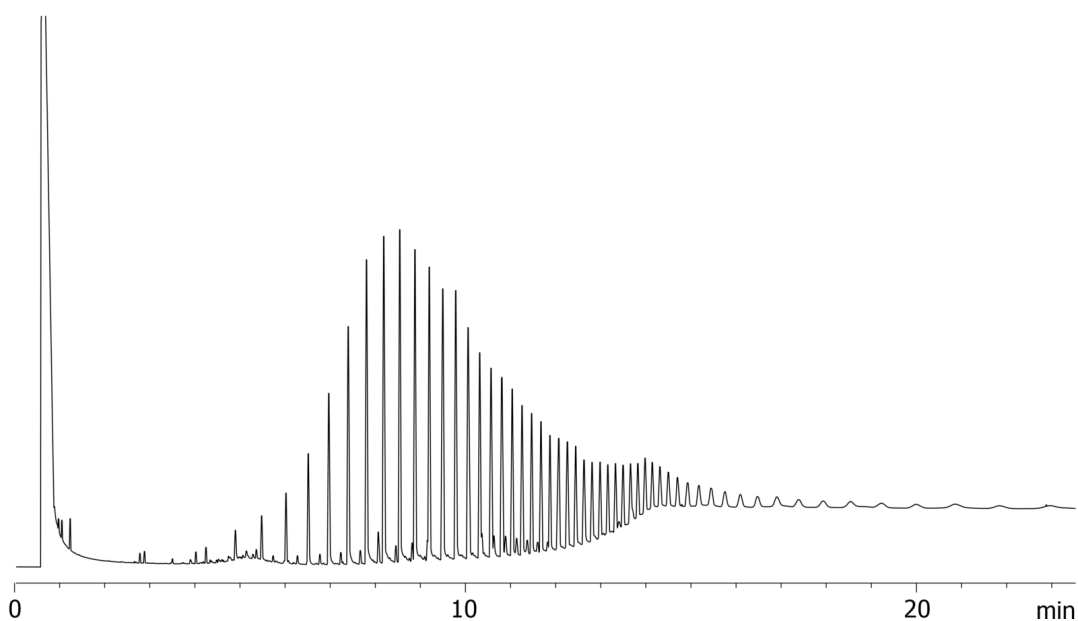


Siloxanes on ZB-5HT**Column:** Zebron[™] ZB-5HT, GC Cap. Column 15 m x 0.32 mm x 0.10 µm, Ea**Phase:** 5% Phenyl 95% Dimethylpolysiloxane**Dimensions:** 15 meters x 0.32 mm x 0.1 µm**Order No:** 7EM-G015-02**Oven Profile:** 60°C to 400°C at 25°C/min for 10min**Carrier Gas:** Constant Flow Helium, 2.9 mL/min**Injection:** On-Column :1 0.1 µL @ 63°C**Detection:** Flame Ionization (FID) (425°C)**Analyst Note:** Extracted septum with carbon disulfide

16122

**Products used in this application:****ANALYTES:****1** Siloxanes