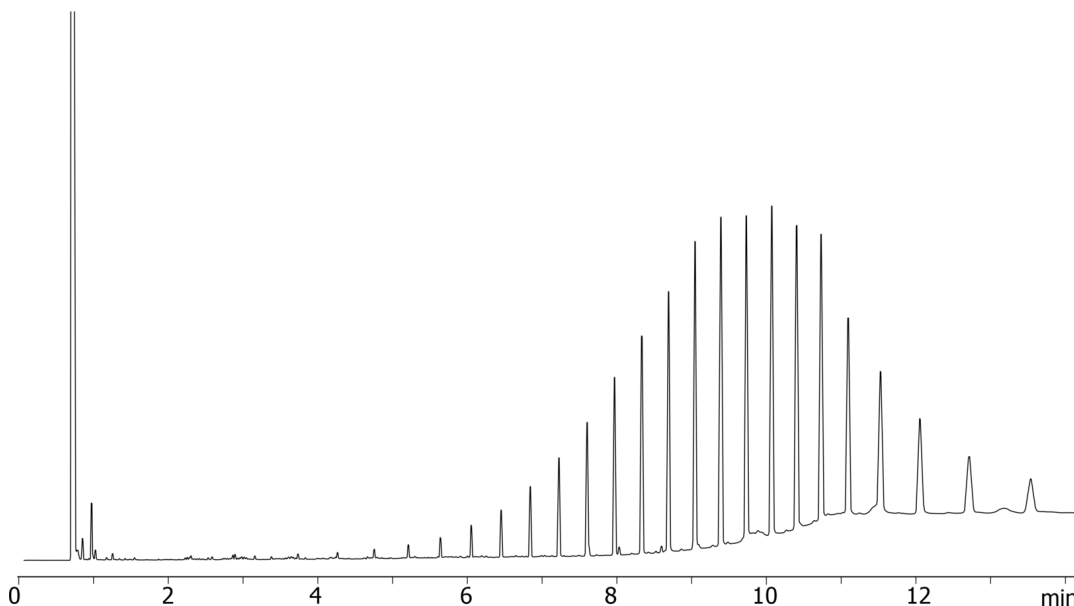


Polywax® 655 on ZB-5HT**Column:** Zebron™ ZB-5HT, GC Cap. Column 15 m x 0.32 mm x 0.10 µm, Ea**Phase:** 5% Phenyl 95% Dimethylpolysiloxane**Dimensions:** 15 meters x 0.32 mm x 0.1 µm**Order No:** 7EM-G015-02**Oven Profile:** 60°C to 380°C at 25°C/min for 5min**Carrier Gas:** Constant Flow Helium, 2.9 mL/min**Injection:** On-Column :1 0.1 µL @ 63°C**Detection:** Flame Ionization (FID) (400°C)**Analyst Note:** Sample was 1.0mg/mL in carbon disulfide. Polywax is a trademark of Petrolite Specialty Polymers Group.

16125

**Products used in this application:****ANALYTES:****1** Polywax® 655