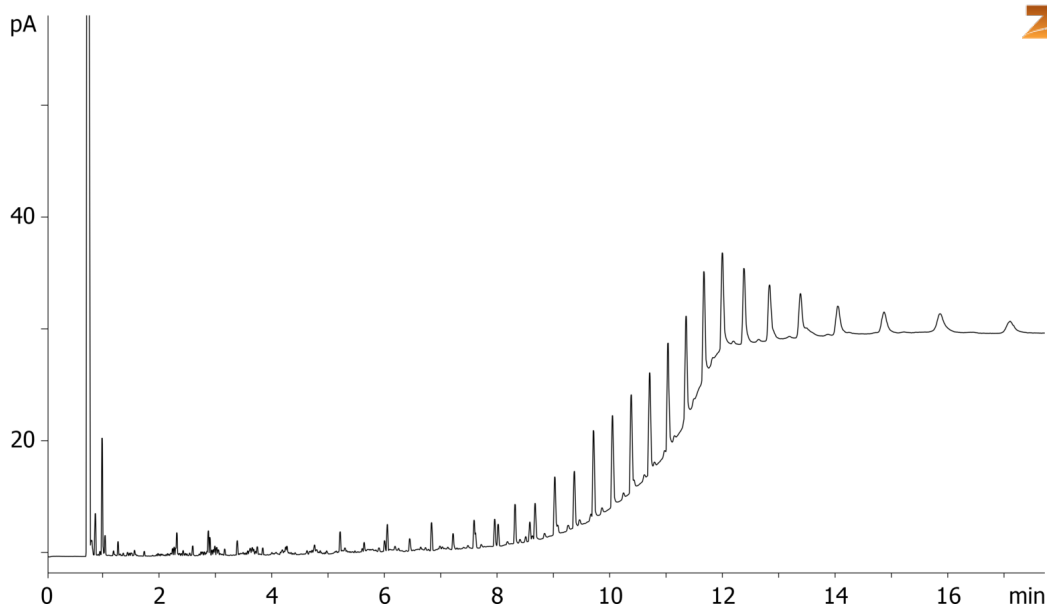


**Polywax® 1000 on ZB-5HT****Column:** Zebtron™ ZB-5HT, GC Cap. Column 15 m x 0.32 mm x 0.10 µm, Ea**Phase:** 5% Phenyl 95% Dimethylpolysiloxane**Dimensions:** 15 meters x 0.32 mm x 0.1 µm**Order No:** 7EM-G015-02**Oven Profile:** 60°C to 400°C at 35°C/min for 8min**Carrier Gas:** Constant Flow Helium, 2.9 mL/min**Injection:** On-Column :1 0.1 µL @ 63°C**Detection:** Flame Ionization (FID) (400°C)**Analyst Note:** Sample was 0.4mg/mL in carbon disulfide. Polywax is a trademark of Petrolite Specialty Polymers Group.**Products used in this application:**

16127

**ANALYTES:****1** Polywax® 1000