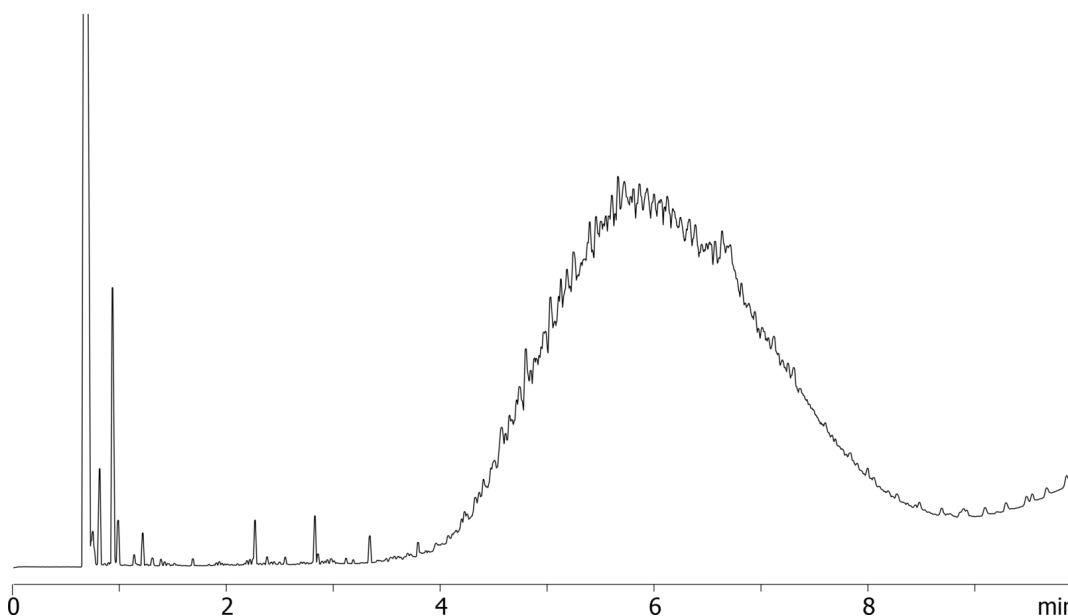


New Motor Oil (10-W30) on ZB-5HT**Column:** Zebron[™] ZB-5HT, GC Cap. Column 15 m x 0.32 mm x 0.10 μ m, Ea**Phase:** 5% Phenyl 95% Dimethylpolysiloxane**Dimensions:** 15 meters x 0.32 mm x 0.1 μ m**Order No:** 7EM-G015-02**Oven Profile:** 60°C to 400°C at 35°C/min for 8 min**Carrier Gas:** Constant Flow Helium, 2.9 mL/min**Injection:** On-Column :1 0.1 μ L @ 63°C**Detection:** Flame Ionization (FID) (400°C)**Analyst Note:** Sample was 1mg/mL in dichloromethane**Products used in this application:**

16128

**ANALYTES:****1** New Motor Oil