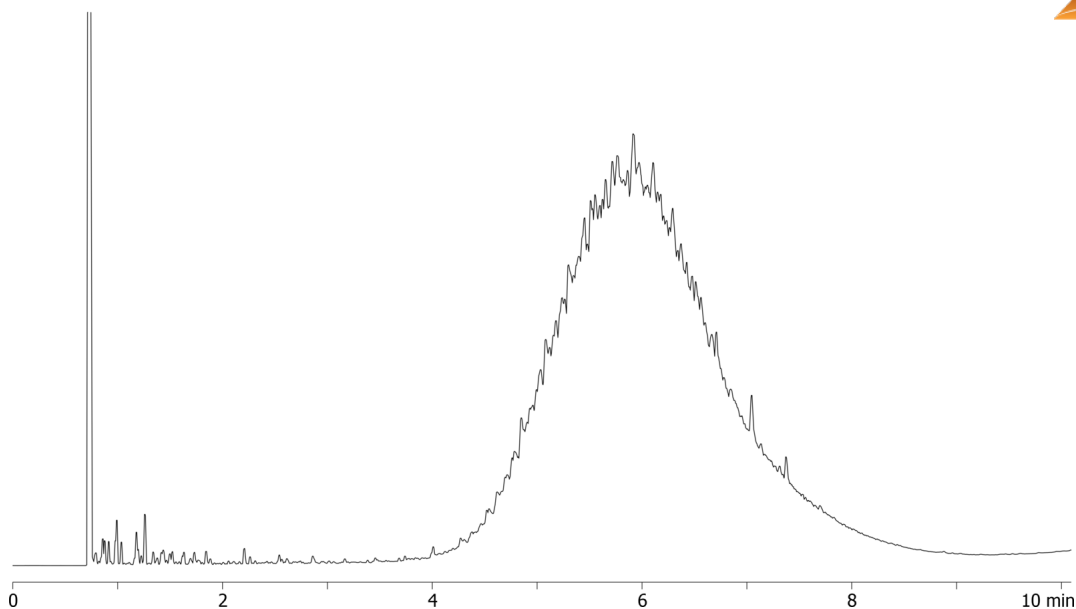


**Used Motor Oil (10-W30) on ZB-5HT****Column:** Zebron<sup>™</sup> ZB-5HT, GC Cap. Column 15 m x 0.32 mm x 0.10  $\mu$ m, Ea**Phase:** 5% Phenyl 95% Dimethylpolysiloxane**Dimensions:** 15 meters x 0.32 mm x 0.1  $\mu$ m**Order No:** 7EM-G015-02**Oven Profile:** 60°C to 400°C at 35°C/min for 8 min**Carrier Gas:** Constant Flow Helium, 2.9 mL/min**Injection:** On-Column :1 0.1  $\mu$ L @ 63°C**Detection:** Flame Ionization (FID) (400°C)**Analyst Note:** Sample was 5mg/mL in dichloromethane**Products used in this application:**

16129

**ANALYTES:****1** Motor oil