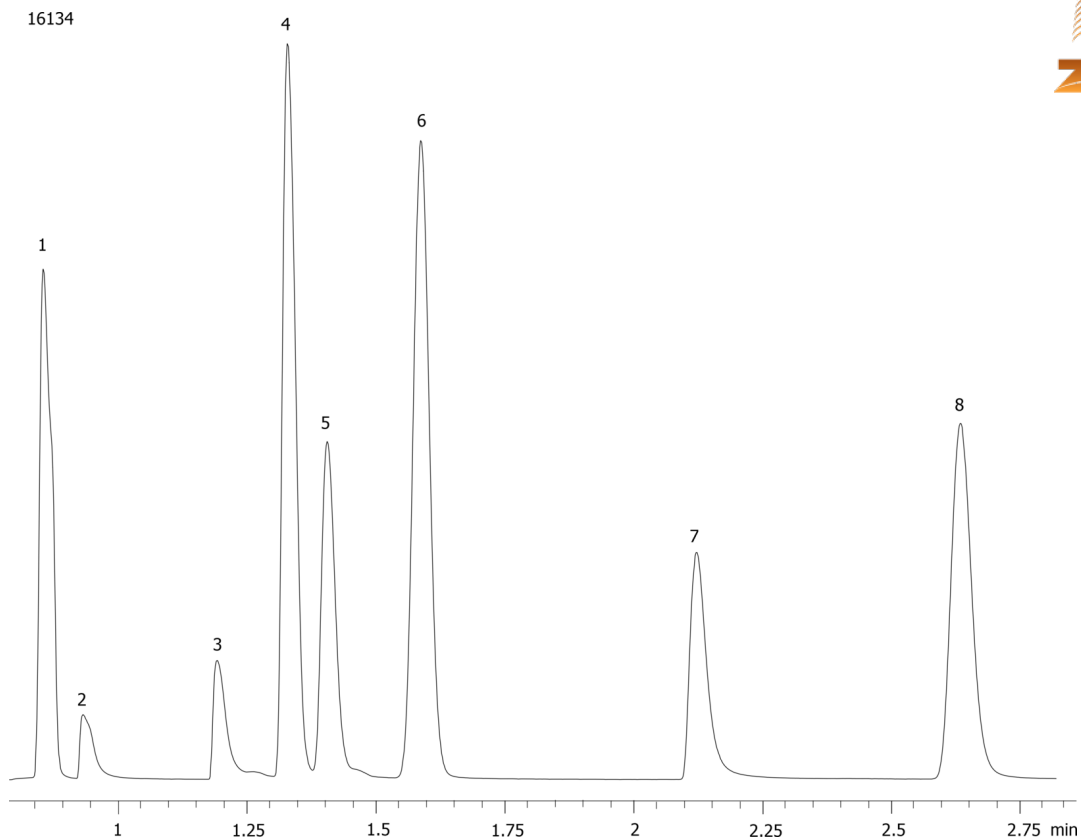


Blood Alcohols on ZB-1701

Column: Zebtron™ ZB-1701, GC Cap. Column 30 m x 0.32 mm x 1.00 µm, Ea
Phase: 14% Cyanopropylphenyl 86% Dimethylpolysiloxane
Dimensions: 30 meters x 0.32 mm x 1 µm
Order No: 7HM-G006-22
Oven Profile: 40°C Isothermal

Carrier Gas: Constant Flow Helium, 7 mL/min
Injection: Split 13:1 200 µL @ 200°C
Detection: Flame Ionization (FID) (220°C)
Analyst Note: Static Headspace: 1mL 50ppm sample in 20mL vial at 70°C for 15 minutes.

**Products used in this application:****ANALYTES:**

- 1 Acetaldehyde
- 2 Methanol
- 3 Ethanol
- 4 Acetone
- 5 Isopropanol
- 6 t-Butanol
- 7 n-Propanol
- 8 2-Butanol

