


European Red List Volatiles on ZB-624**Column:** Zebtron[™] ZB-624, GC Cap. Column 20 m x 0.18 mm x 1.00 µm, Ea**Phase:** 6% Cyanopropylphenyl 94% Dimethylpolysiloxane**Dimensions:** 20 meters x 0.18 mm x 1 µm**Order No:** 7FD-G005-22**Oven Profile:** 90°C Isothermal**Carrier Gas:** Constant Flow Helium, 1.5 mL/min**Injection:** Split 50:1 1 µL @ 250°C**Detection:** Flame Ionization (FID) (260°C)**Analyst Note:** Static Headspace: 8µL of 12.5% volatiles mix sample in a 0.18mL limited volume insert inside a 2.0mL vial.**Products used in this application:****ANALYTES:**

- 1** 1,1-Dichloroethylene
- 2** 1,1-Dichloroethane
- 3** Chloroform
- 4** 1,1,1-Trichloroethane
- 5** Carbon tetrachloride
- 6** 1,2-Dichloroethane
- 7** Trichloroethylene
- 8** Tetrachloroethylene

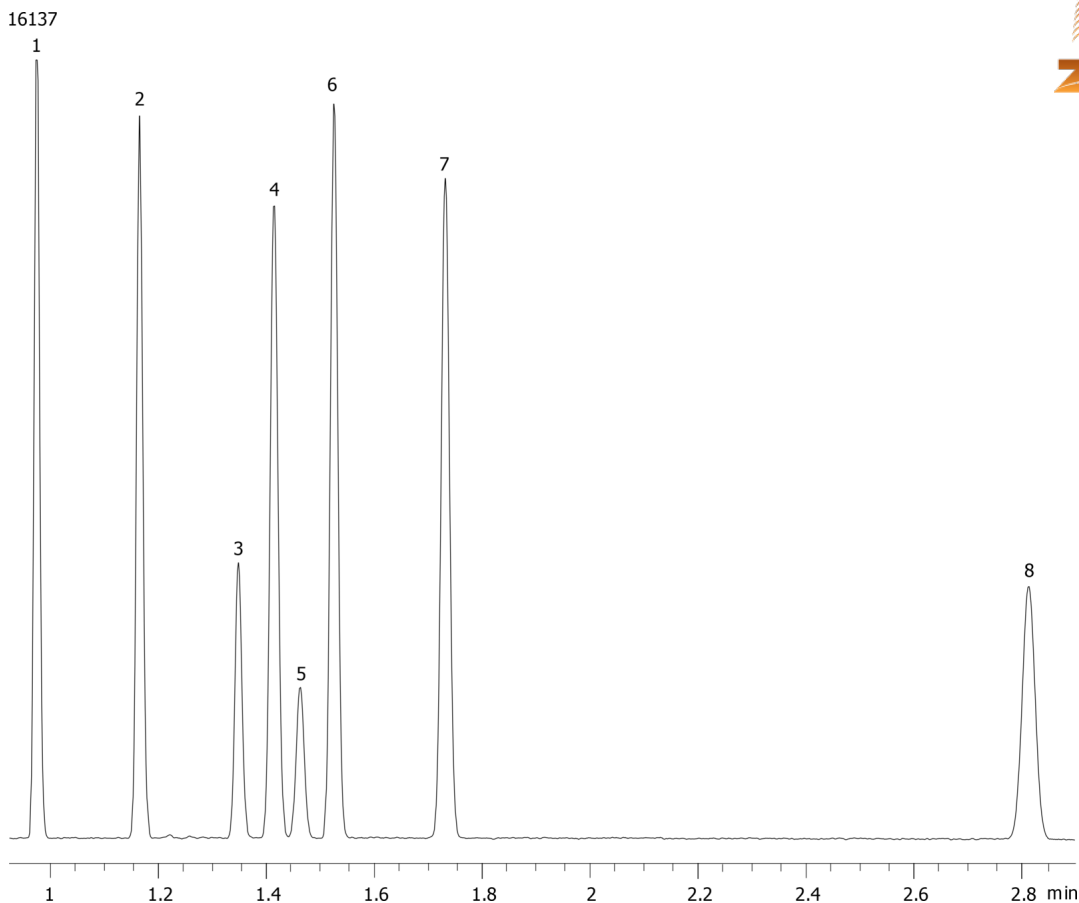