


Low Breakdown of DDT on ZB-MR1**Column:** Zebron[™] ZB-MultiResidue[™]-1, GC Cap. Column 30 m x 0.32 mm x 0.50 μ m, Ea**Phase:****Dimensions:** 30 meters x 0.32 mm x 0.5 μ m**Order No:** 7HM-G016-17**Oven Profile:** 100 °C for 0.5 min to 340 °C at 15 °C/min for 2 min**Carrier Gas:** Constant Flow Helium, 1.6 mL/min**Injection:** Splitless :1 1 μ L @ 250°C**Detection:** Electron Capture (ECD) (350°C)**Analyst Note:** Analyte was 250 pg on column.

16165

**Products used in this application:****ANALYTES:**

- 1** DDT
- 2** DDD
- 3** DDE

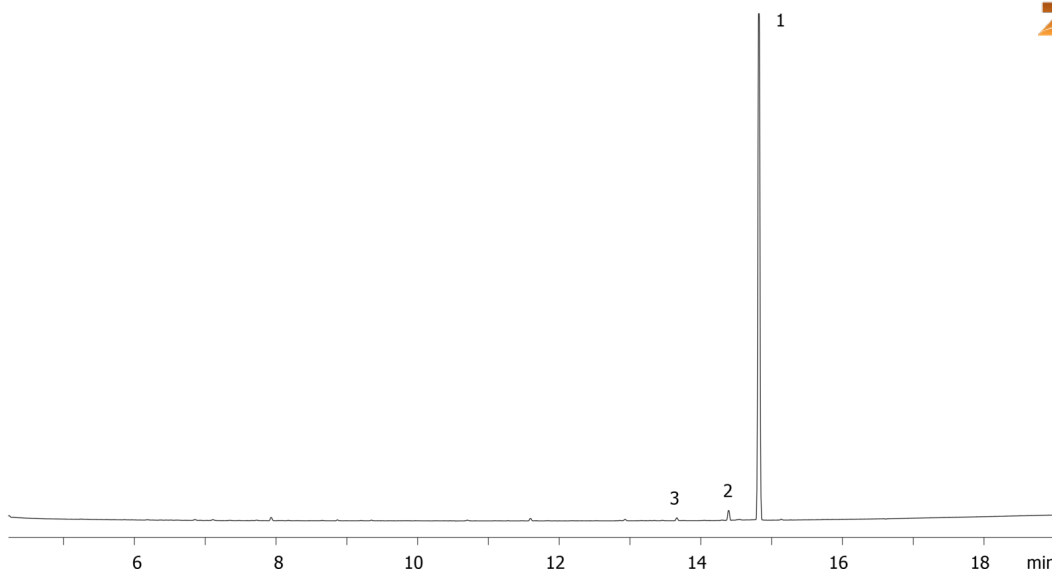