


Aroclor 1242 on ZB-MR1**Column:** Zebron[™] ZB-MultiResidue[™]-1, GC Cap. Column 30 m x 0.32 mm x 0.50 μ m, Ea**Phase:****Dimensions:** 30 meters x 0.32 mm x 0.5 μ m**Order No:** 7HM-G016-17**Oven Profile:** 120 °C hold 1min to 300 °C at 9 °C/min hold 10 min**Carrier Gas:** Constant Pressure Hydrogen, 2.3 mL/min**Injection:** Split 50:1 1 μ L @ 210°C**Detection:** Electron Capture (ECD) (310°C)

16167

**Products used in this application:****ANALYTES:**

- 1** Tetrachloro-m-xylene
- 2** Aroclor 1242
- 3** Decachlorobiphenyl

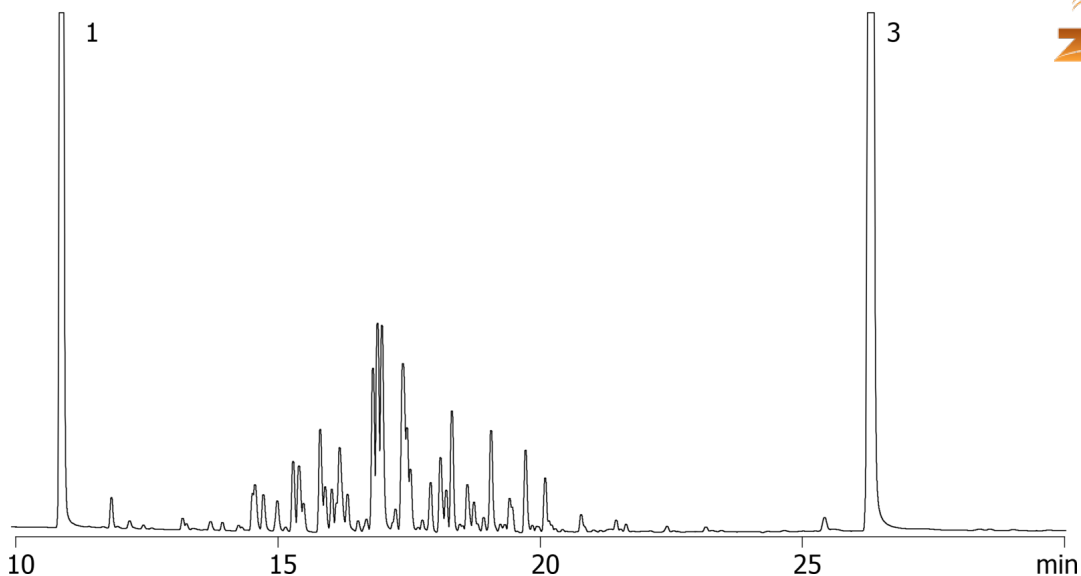