


Aroclor 1268 on ZB-MR2**Column:** Zebtron[™] ZB-MultiResidue[™]-2, GC Cap. Column 30 m x 0.32 mm x 0.25 μ m, Ea**Phase:****Dimensions:** 30 meters x 0.32 mm x 0.25 μ m**Order No:** 7HM-G017-11**Oven Profile:** 120 °C hold 1min to 300 °C at 9 °C/min hold 10 min**Carrier Gas:** Constant Pressure Hydrogen, 2.3 mL/min**Injection:** Split 50:1 1 μ L @ 210°C**Detection:** Electron Capture (ECD) (310°C)

16174

**Products used in this application:****ANALYTES:**

- 1 Tetrachloro-m-xylene
- 2 Aroclor 1268
- 3 Decachlorobiphenyl

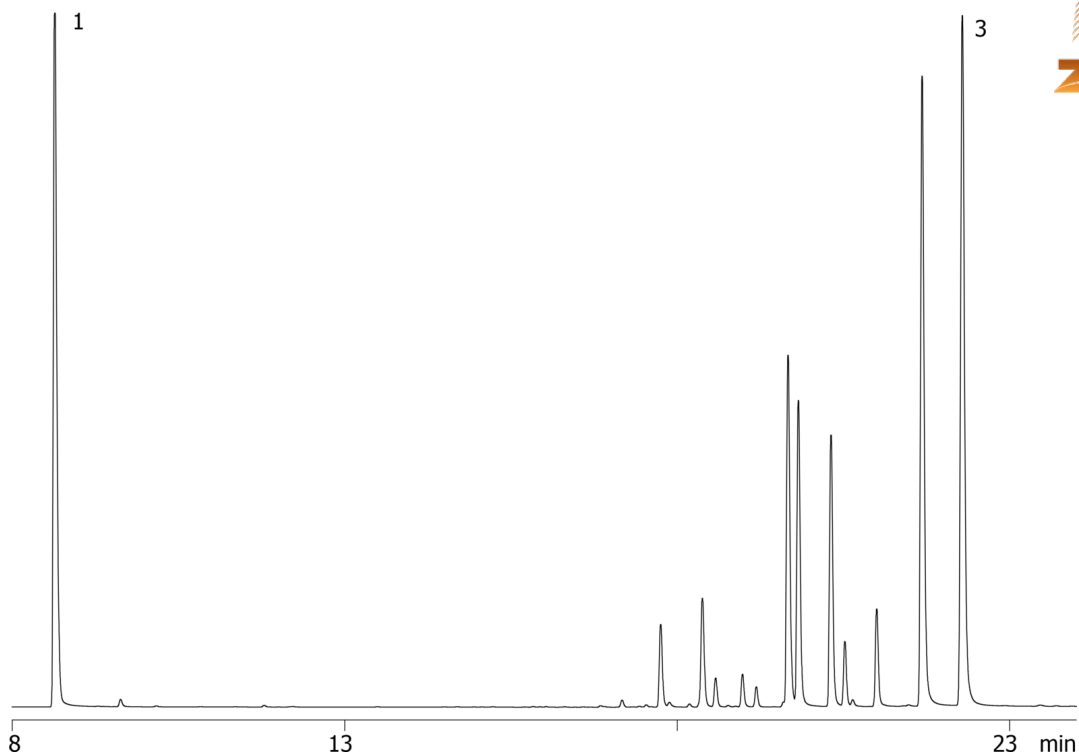