

GC Application

ID No.: 16178

Aroclor 1242: DIN 51527 on ZB-MR2

Column: Zebron[™] ZB-MultiResidue[™]-2, GC Cap. Column 30 m x 0.32 mm x 0.25 μm, Ea

Phase:

Dimensions: 30 meters x 0.32 mm x 0.25 μm

Order No: 7HM-G017-11


Oven Profile: 100°C for 0.5min to 200°C at 25°C/min to 320°C at 6°C/min for 2min

Carrier Gas: Constant Flow Helium, 1.5 mL/min

Injection: Splitless :1 1 μL @ 250°C

Detection: Electron Capture (ECD) (350°C)

16178

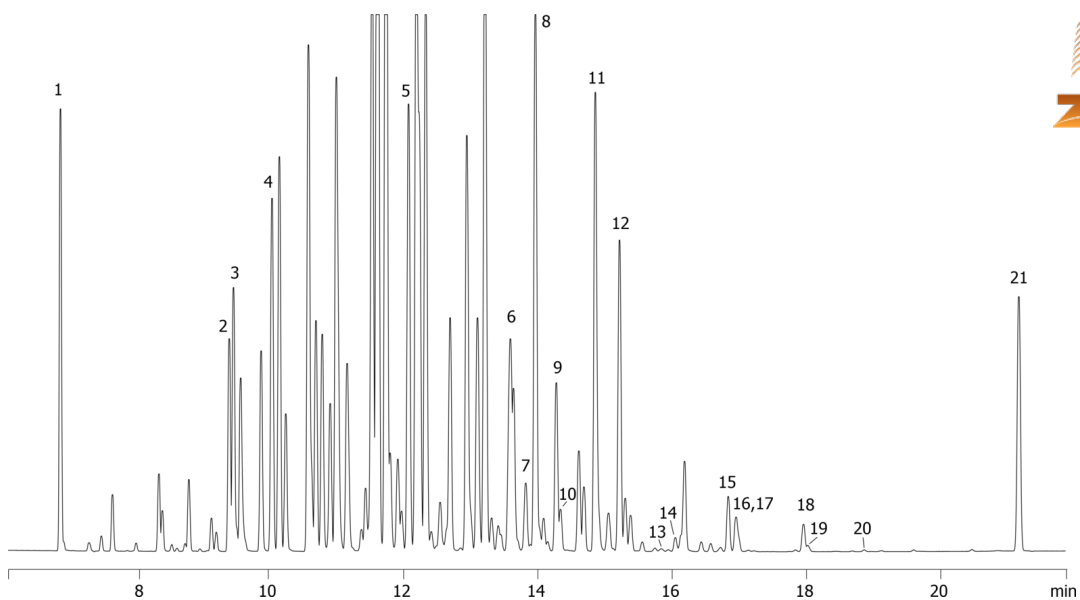

Products used in this application:



Aroclor 1242: DIN 51527 on ZB-MR2

ANALYTES:

- 1 TCMX
- 2 BZ# 31
- 3 BZ# 28
- 4 BZ# 52
- 5 BZ# 101
- 6 BZ# 77
- 7 BZ# 123
- 8 BZ# 118
- 9 BZ# 153
- 10 BZ# 114
- 11 BZ# 105
- 12 BZ# 138
- 13 BZ# 126
- 14 BZ# 167
- 15 BZ# 156
- 16 BZ# 180
- 17 BZ# 157
- 18 BZ# 170
- 19 BZ# 169
- 20 BZ# 189
- 21 DCB

