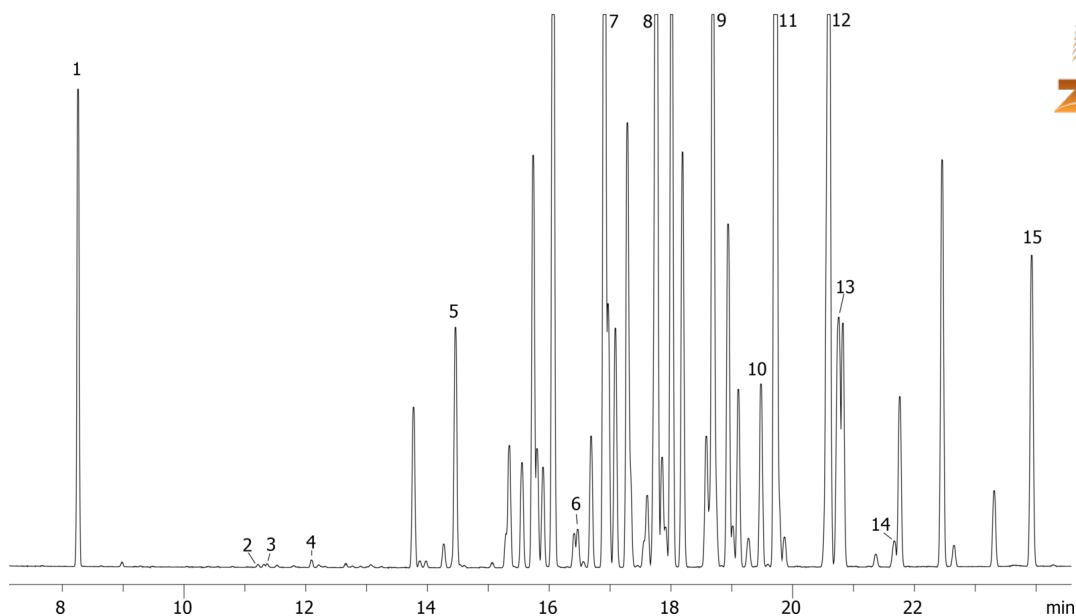


Aroclor 1260: DIN 51527 on ZB-MR1**Column:** Zebron[™] ZB-MultiResidue[™]-1, GC Cap. Column 30 m x 0.32 mm x 0.50 μ m, Ea**Phase:****Dimensions:** 30 meters x 0.32 mm x 0.5 μ m**Order No:** 7HM-G016-17**Oven Profile:** 100 °C for 0.5 min to 200 °C at 25 °C/min to 320 °C at 6 °C/min for 2 min**Carrier Gas:** Constant Flow Helium, 1.5 mL/min**Injection:** Splitless :1 1 μ L @ 250°C**Detection:** Electron Capture (ECD) (350°C)

16179

**Products used in this application:**

Aroclor 1260: DIN 51527 on ZB-MR1

ANALYTES:

- 1 TCMX
- 2 BZ# 31
- 3 BZ# 28
- 4 BZ# 52
- 5 BZ# 101
- 6 BZ# 118
- 7 BZ# 153
- 8 BZ# 138
- 9 BZ# 167
- 10 BZ# 156
- 11 BZ# 180
- 12 BZ# 170
- 13 BZ# 169
- 14 BZ# 189
- 15 DCB

